(12) INTERNATIONAL APPLICATION PUBLISHED UNDER THE PATENT COOPERATION TREATY (PCT)

(19) World Intellectual Property Organization

International Bureau

(43) International Publication Date 11 August 2011 (11.08.2011)





(10) International Publication Number WO 2011/097503 A2

(51) International Patent Classification: *C12Q 1/68* (2006.01)

(21) International Application Number:

PCT/US2011/023777

(22) International Filing Date:

4 February 2011 (04.02.2011)

(25) Filing Language:

English

(26) Publication Language:

English

(30) Priority Data:

61/302,059 5 February 2010 (05.02.2010)

US

(71) Applicant (for all designated States except US): QUEST DIAGNOSTICS INVESTMENTS INCORPORATED [US/US]; 300 Delaware Avenue, Wilmington, DE19899 (US).

(72) Inventors; and

(71) Applicants: HANTASH, Feras [CA/US]; 25071 Bellota, Mission Viejo, CA 92692 (US). SUN, Weimin [US/US]; 23 Cambridge, Irvine, CA 92620 (US). TSAO, David, C. [US/US]; 1720 Orchard Hill Lane, Hacienda Heights, CA 91745 (US). ROOT, Dana, Marie [US/US];

8902 Heil Ave., Apt. #3, Westminster, CA 92683 (US). STROM, Charles, M. [US/US]; 2939 Calle Gaucho, San Clemente, CA 92673 (US).


(74) Agent: KUHLMANN, Anthony, C.; Foley & Lardner LLP, P.O. Box 80278, San Diego, CA 92138 (US).


(81) Designated States (unless otherwise indicated, for every kind of national protection available): AE, AG, AL, AM, AO, AT, AU, AZ, BA, BB, BG, BH, BR, BW, BY, BZ, CA, CH, CL, CN, CO, CR, CU, CZ, DE, DK, DM, DO, DZ, EC, EE, EG, ES, FI, GB, GD, GE, GH, GM, GT, HN, HR, HU, ID, IL, IN, IS, JP, KE, KG, KM, KN, KP, KR, KZ, LA, LC, LK, LR, LS, LT, LU, LY, MA, MD, ME, MG, MK, MN, MW, MX, MY, MZ, NA, NG, NI, NO, NZ, OM, PE, PG, PH, PL, PT, RO, RS, RU, SC, SD, SE, SG, SK, SL, SM, ST, SV, SY, TH, TJ, TM, TN, TR, TT, TZ, UA, UG, US, UZ, VC, VN, ZA, ZM, ZW.

(84) Designated States (unless otherwise indicated, for every kind of regional protection available): ARIPO (BW, GH, GM, KE, LR, LS, MW, MZ, NA, SD, SL, SZ, TZ, UG, ZM, ZW), Eurasian (AM, AZ, BY, KG, KZ, MD, RU, TJ, TM), European (AL, AT, BE, BG, CH, CY, CZ, DE, DK, EE, ES, FI, FR, GB, GR, HR, HU, IE, IS, IT, LT, LU, LV, MC, MK, MT, NL, NO, PL, PT, RO, RS, SE, SI, SK,

[Continued on next page]

(54) Title: METHOD TO DETECT REPEAT SEQUENCE MOTIFS IN NUCLEIC ACID

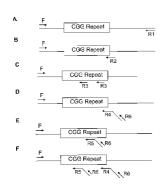

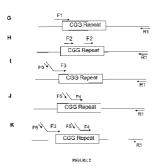

(57) Abstract: Methods for determining the presence or absence of expansion of CGG repeat sequence in the FMR1 gene presence or absence of expansion of CCG repeat sequence in the FMR2 gene are provided. The methods are useful in identifying an individual with normal/intermediate, versus premutation or full mutation allele of FMR1 gene and FMR2 gene due to the expansion of CGG repeats and CCG repeats in the 5'-untranslated region respectively. The methods are also useful for screening newborns for fragile X syndrome or for screening women to determine heterozygosity status with full premutation of the CCG repeat tract. The methods are also useful in estimating the premutation and full mutation carrier frequency and estimating the prevalence of FXTAS AND FXPOI in a population. The methods are simple, rapid and require small amount of sample.



SM, TR), OAPI (BF, BJ, CF, CG, CI, CM, GA, GN, GQ, **Published**: GW, ML, MR, NE, SN, TD, TG).

without international search report and to be republished upon receipt of that report (Rule 48.2(g))

METHOD TO DETECT REPEAT SEQUENCE MOTIFS IN NUCLEIC ACID

CROSS-REFERENCE TO RELATED PATENT APPLICATIONS

[0001] This application claims the priority of U.S. Provisional Application No. 61/302,059 filed February 5, 2010 which is incorporated herein by reference in its entirety.

FIELD OF THE INVENTION

[0002] This invention relates to methods of detecting genetic mutations characterized by an expansion of tandem repeats such as CGG repeats in the *FMR1* gene and CCG repeats in FMR2 gene.

BACKGROUND

[0003] A tandem repeat in DNA represents two or more contiguous approximate copies of a pattern of nucleotides. Tandem repeats have been shown to be associated with a variety of human diseases. Dramatic expansion of trinucleotide repeats has been associated with such diseases as fragile-X mental retardation (see Verkerk, et al., (1991) Cell, 65, 905-914), Huntington's disease (see Huntington's Disease Collaborative Research Group. (1993) Cell, 72, 971-983), myotonic dystrophy (see Fu, et al., (1992) Science, 255, 1256-1258), spinal and bulbar muscular atrophy (see La Spada, et al., (1991) Nature, 352, 77-79) and Friedreich's ataxia (see Campuzano, et al., (1996) Science, 271, 1423-1427).

[0004] Fragile X syndrome (FXS, OMIM # 300624) is the most common inherited mental retardation syndrome, affecting ~1:4000 males and ~1:8000 females. The disease is caused by the expansion of the trinucleotide CGG repeat in the 5'-untranslated (UTR) region of the fragile X mental retardation protein 1 (*FMR1*) gene. Methylation of the CGG tract leads to silencing of expression of the *FMR1* gene. The American College of Medical Genetics (Maddalena et al. Genet Med 2001; 3(3):200-5; Sherman et al. Genet Med 2005;7(8):584-7) defines a normal repeat length as between 5-44. Intermediate alleles of between 45-54 repeats show almost no expansion. Premutation alleles have 55-200 CGG, while full mutation alleles have >200 repeats. As an X-linked dominant disease, FXS is diagnosed more frequently in males than females, but females can present with milder symptoms, or with FXS due to

skewed X-inactivation. Thirty-forty percent of male carriers of a premutation allele will suffer from Fragile X-associated Tremor and Ataxia (FXTAS) by age 50 (Hagerman et al. Ment Retard Dev Disabil Res Rev 2004; 10(1): 25-30; Hagerman et al. Am J Hum Genet 2004; 74(5): 1051-6) Various studies showed that the larger the repeat size the stronger the chance of presenting with FXTAS, but the penetrance i.e. the chance of being affected by FXTAS seems to be around 33% for male carriers of alleles ~70 repeats or higher. Some female premutation carriers are affected with premature ovarian insufficiency (FXPOI; also know as premature ovarian failure) (Allingham-Hawkins et al. Am J Med Genet 1999;83(4):322-5; Hundscheid et al. Hum. Reprod. 2001; 13(3): 457-462). Female premutation carriers showed overall higher levels of follicular stimulating hormone than normal or full mutation carriers and decreased levels of anti-Müllerian hormone, Inhibin A and inhibin B, all indicators of ovarian decline (Hundscheid et al. Hum. Reprod. 2001; 13(3): 457-462). The reported penetrance i.e. the chance of being affected by FXPOI in female premutation carriers is 20-28% of larger alleles.

[0005] Several studies have examined the carrier frequency of premutation and/or full mutation alleles in FMR1 gene in different populations, and as in any other genetic disease, the frequency varies by population and region; 1:113 in Israel (Toledano-Alhadef et al. Am J Hum Genet 2001; 69(2):351-60), 1:259 in Quebec (Rousseau et al. Am J Hum Genet 1995;57(5):1006-18), in Finland the frequency was 1:246, and in the U.S. state of Georgia it was 1:436. The frequency was found to be between 1:257-1:382 in the U.S. population in general (Cronister et al. Genet Med 2005;7(4):246-50). Methods of amplifying the CGG repeat tract in the FMR1 gene has previously been reported. See, Fu et al. Cell 1991; 67: 1047-1058; Erster et al. Hum Genet. 1992; 90: 55-61; Pergolizzi et al. Lancet. 1992; 339: 271-272; Warner et al. J Med Genet. 1996; 33: 1022-1026; Daniels et al. J Assist Reprod Genet. 1996; 13(2): 163-169; Strelnikov et al. Hum Mutat. 1999; 13(2): 166-169; Zhou et al. J Med Genet. 2004; 41(4): e45; Saluto et al. J Mol Diagn. 2005; 7(5): 605-612; Zhou et al. Clin Chem. 2006; 52(8): 1492-1500; Dahl et al. Clin Chem. 2007; 53(4): 790-793; Strom et al. Genet Med. 2007; 9(4): 199-207; Khaniani et al. Mol Cytogenet. Apr 8; 1:5; Tassone et al. J Mol Diagn. 2008; 10(1): 43-49; Filipovic-Sadic et al. Clin Chem. January 7 2010 (epub); Lyon et al. J Mol. Diagn. 2010; 12(4): 505-511; Hantash et al. Genet Med. Feb 17 2010 (epub); US Patent 5,658,764, 6,200,747, and 6,114,150, and US Patent Application Publication 2008/0113355.

[0006] Fragile X E (FRAXE) mental retardation (OMIM 309548) is associated to a fragile site localized in chromosome Xq28 and is the cause of a non-syndromic X-linked mental retardation affecting 1/50,000 newborn males. The disorder is due to the silencing of the Fragile X Mental Retardation 2 (FMR2) gene, as a consequence of a CCG expansion located upstream to this gene. FRAXE alleles can be divided into four categories: normal (6-30 CCG repeats), intermediate (31-60 CCG repeats), premutation (61-200 CCG repeats), and full mutation (over 200 CCG repeats). See, Gecz et al. Nat. Genet. 1996;13: 105–108; Gu et al. Nat. Genet. 1996;13:109–113; Knight et al. Cell. 1993;74:127–134; and Strelnikov et al. Hum Mutat. 1999;13(2):166-9.

SUMMARY OF THE INVENTION

[0007] Described is a nucleic acid amplification method for detecting whether a gene associated with a tandem triplet repeat tract has a normal, or an expanded triplet repeat track. The expanded triplet repeat tract may include a premutation or may include a full mutation with a greater number of repeats than the premutation. In one embodiment, amplification is by the polymerase chain reaction using a pair of primers. The amplified triplet repeat tract and any portions thereof are separated according to size. In one embodiment, the method of separation is electrophoresis, which preferably is capillary gel electrophoresis.

[0008] In one aspect, the method is used to determine the presence or absence of an expansion of the CGG repeat tract in the 5'-untranslated region (5'-UTR) of fragile X related mental retardation gene 1 (FMR1) in an individual. The method involves amplification using at least a primer pair wherein one or both primers of the primer pair hybridize to and prime in the CGG repeat tract so as to generate amplicons that contain all or portions of the CGG repeat segment which contains less than four consecutive CCG triplet repeats and hybridizes to at least a portion of said CGG repeat tract. Alternatively or in addition, the primer pair can include an upstream primer that comprises a CGG repeat segment which contains less than four consecutive CGG triplet repeats and hybridizes to at least a portion of the complement of said CGG repeat tract. The amplicons generated are then separated according to size, wherein the size(s) of the amplicons indicates the presence or absence of expansion of the CGG repeat tract in said 5'-UTR of the FMR1 gene.

[0009] In some embodiments, when the primer pair for amplifying CGG repeat tract in the 5'-UTR of FMR1 gene includes a downstream primer that comprises less than four consecutive CCG triplet repeats and hybridizes to at least a portion of the CGG repeat tract, the upstream primer is complementary to a region upstream of the CGG repeat tract. Likewise, when the primer pair includes an upstream primer that comprises less than four consecutive CGG triplet repeats and hybridizes to at least a portion of the complement of the CGG repeat tract, the downstream primer is complementary to sequence downstream of the CGG tract.

[0010] In another aspect, the method is used to determine the presence or absence of an expansion of the CCG repeat tract in the 5'-untranslated region (5'-UTR) of fragile X related mental retardation gene 2 (FMR2) in an individual. The method involves amplification using at least a primer pair wherein one or both primers of the primer pair hybridize to and prime in the CCG repeat tract so as to generate amplicons that contain all or portions of the CCG repeat tract. The primer pair can include a downstream primer that comprises one or more CGG triplet repeats and hybridizes to at least a portion of said CCG repeat tract. In some embodiments, the downstream primer may comprise two or more consecutive CGG triplet repeats. Alternatively or in addition, the primer pair can include an upstream primer that comprises one or more CCG repeat segment and hybridizes to at least a portion of the complement of said CCG repeat tract. In some embodiments, the upstream primer may comprise two or more consecutive CCG triplet repeats. The amplicons generated are then separated according to size, wherein the size(s) of the amplicons indicates the presence or absence of expansion of the CCG repeat tract in said 5'-UTR of the FMR2 gene.

[0011] In some embodiments, when the primer pair for amplifying CCG repeat tract in the 5'-UTR of FMR2 gene includes a downstream primer that comprises one or more consecutive CGG triplet repeats and hybridizes to at least a portion of the CCG repeat tract, the upstream primer is complementary to a region upstream of the CCG repeat tract. Likewise, when the primer pair includes an upstream primer that comprises one or more consecutive CCG triplet repeats and hybridizes to at least a portion of the complement of the CCG repeat tract, the downstream primer is complementary to sequence downstream of the CCG tract.

[0012] The primers can include 5' tails that can function as a priming site and have sequence with little to no sequence homology to the FMR1 or FMR2 gene. Thus, an additional primer (e.g. M13 sequence) can be used which primes in the tail segment of one or both primers of the initial primer to boost the amount of amplicons generated.

[0013] The CGG or CCG repeat segments of the primers can have various combinations of triplets such as are shown in FIG. 15 for primers comprising CCG triplets and FIG. 16 for primers comprising CGG triplets. In one embodiment, the primer that has the CGG or CCG repeat segment also includes sequence that flanks the CGG or CCG repeat track so as to allow the primer to hybridize at the junction between the CGG repeat tract or CCG repeat tract and flanking sequence. For example, in the case of FMR1 gene, the "junctional" downstream primer has CCG repeat sequence at the 3' end of the primer which sequence includes sequence upstream from the 3' end of the repeat tract and continues into the flanking sequence that is directly downstream of the repeat tract. Likewise, the "junctional" form of the upstream primer has at its 3' end CGG repeat sequence that includes sequence downstream from the 5' end of the repeat tract and continues upstream into the flanking sequence located 5' to the repeat tract. As already discussed, these junctional primers can include a tail for priming by an additional primer.

[0014] In some embodiments, for amplifying CGG repeat tract in the 5'-UTR of FMR1 gene the downstream primer comprises or consists essentially of a sequence selected from the group consisting of SEQ ID NOs: 3-6, SEQ ID NOs 14-40, SEQ ID NOs 71-74 and the upstream primer comprises or consists essentially of a sequence selected from the group consisting of SEQ ID NOs: 41-69. In some embodiments, for amplifying CGG repeat tract in the 5'-UTR of FMR1 gene the CGG or CCG repeat segments of the primer can have less than 3 consecutive CGG or CCG repeats. In some embodiments, the downstream primer CCG hybridizing sequence has two sets of two consecutive CCG triplets (separated by at least one base) and the upstream primer has as the CGG hybridizing sequence two sets of two consecutive CGG triplets (separated by at least one base). In these embodiments, the downstream primer comprises or consists essentially of a sequence selected from the group consisting of SEQ ID NOs: 15-22 and 29-36 and the upstream primer comprises or consists essentially of a sequence selected from the group consisting of SEQ ID NOs: 42-49 and 56-66.

gene the upstream primer comprises or consists essentially of a sequence selected from the group consisting of SEQ ID NOs: 14-40, 71-74, and 76 and the downstream primer comprises or consists essentially of a sequence selected from the group consisting of SEQ ID NOs: 41-69 and 77. In some embodiments, for amplifying CCG repeat tract in the 5'-UTR of FMR2 gene, the CGG or CCG repeat segments of the primer can have less than 3 consecutive CGG or CCG repeats. In some embodiments, the upstream primer has two sets of two consecutive CCG triplets (separated by at least one base) and the downstream primer has two sets of two consecutive CGG triplets (separated by at least one base). In these embodiments, the upstream primer comprises or consists essentially of a sequence selected from the group consisting of SEQ ID NOs: 15-22 and 29-36 and the downstream primer comprises or consists essentially of a sequence selected from the group consists essentially of a sequence selected from the group consists essentially of a sequence selected from the group consists essentially of a sequence selected from the group consists essentially of a sequence selected from the group consists essentially of a sequence selected from the group consists essentially of a sequence selected from the group consists essentially of a sequence selected from the group consists essentially of a sequence selected from the group consists essentially of a sequence selected from the group consists essentially of a sequence selected from the group consists essentially of a sequence selected from the group consists essentially of a sequence selected from the group consists essentially of a sequence selected from the group consists essentially of a sequence selected from the group consists essentially of a sequence selected from the group consists essentially of a sequence selected from the group consists essentially of a sequence selected from the group consists essentially essentially essentially essentially essentially essentially essentially essentially

[0016] In accordance with the method, amplification of extended triplet repeats including premutation and full mutation following size separation of the amplicons yields a tapering or stutter pattern of amplicon sizes. In some embodiments, in the case of a full mutation in FMR1 gene, a lone peak is visible at about 1050 bp on a capillary gel.

[0017] In some embodiments, the signal for the full mutation lone peak can be enhanced by using a deoxy GTP analogue such as 7-deaza-2'-deoxy GTP in the amplification reaction mixture. In some embodiments, the dGTP analogue substitutes entirely for dGTP in the reaction mixture.

[0018] In some embodiments, the amplification is cyclic and the ramp rate between cycles is selected at about 60-70% (about 0.4° C to 0.6° C/sec) for the denaturation and/or extension steps and about 20-30% (about 0.4° C to 0.6° C/sec) for the annealing step.

[0019] In another aspect, the method of triplet primed amplification is used to identify a female that is homozygous normal versus a female that is heterozygous and having a normal allele and an allele with an expanded CGG tract length in the 5'-UTR of FMR1 gene. The above described methods may be used in this instance to obtain amplicons which are separated according to size so as to estimate the number of CGG repeats in said amplicons, wherein detection of only a normal number of CGG repeats indicates a homozygous normal female and detection of both amplicons with a normal number of CGG repeats and detection

of amplicons with an abnormally high number of CGG repeats indicates a female heterozygous and having an allele with an expanded CGG tract length. In some embodiments, the normal number of CGG repeats is between 4 to 44 repeats while an abnormally high number of repeats is greater than 54 repeats. In some embodiments, an expanded CGG tract length representing a premutation has CGG repeats of between 55-200 while in other embodiments, an expanded CGG tract length representing a full mutation has greater than 200 CGG triplet repeats.

[0020] In another aspect, the method of triplet primed amplification is used to identify a female that is homozygous normal versus a female that is heterozygous and having a normal allele and an allele with an expanded CCG tract length in the 5'-UTR of FMR2 gene. The above described methods may be used in this instance to obtain amplicons which are separated according to size so as to estimate the number of CCG repeats in said amplicons, wherein detection of only a normal number of CCG repeats indicates a homozygous normal female and detection of both amplicons with a normal number of CCG repeats and detection of amplicons with an abnormally high number of CCG repeats indicates a female heterozygous and having an allele with an expanded CCG tract length. In some embodiments, the normal number of CCG repeats is between 6 to 30 repeats while an abnormally high number of repeats is greater than 60 repeats. In some embodiments, an expanded CCG tract length representing a premutation has CCG repeats of between 61-200 while in other embodiments, an expanded CCG tract length representing a full mutation has greater than 200 CCG triplet repeats.

[0021] In yet a further aspect, the method of triplet primed amplification is used to screen newborns to determine if they have an expanded CGG tract in the FMR1 gene. This method can be performed as for the method above for screening females to determine normal homozygous versus heterozygous with a expanded CGG repeat tract.

[0022] "Primer" as used herein refers to an oligonucleotide that is capable of acting as a point of initiation of synthesis when placed under conditions in which primer extension is initiated (e.g., primer extension associated with an application such as PCR). The primer is complementary to a target nucleotide sequence and it hybridizes to a substantially complementary sequence in the target and leads to addition of nucleotides to the 3′-end of the primer in the presence of a DNA or RNA polymerase. The 3′- nucleotide of the primer should generally be complementary to the target sequence at a corresponding nucleotide

position for optimal expression and amplification. An oligonucleotide "primer" may occur naturally, as in a purified restriction digest or may be produced synthetically. The term "primer" as used herein includes all forms of primers that may be synthesized including peptide nucleic acid primers, locked nucleic acid primers, phosphorothioate modified primers, labeled primers, and the like.

- [0023] Primers are typically between about 10 and about 100 nucleotides in length, preferably between about 15 and about 60 nucleotides in length, more preferably between about 20 and about 50 nucleotides in length, and most preferably between about 25 and about 40 nucleotides in length. In some embodiments, primers can be at least 8, at least 12, at least 16, at least 20, at least 25, at least 30, at least 35, at least 40, at least 45, at least 50, at least 55, at least 60 nucleotides in length. An optimal length for a particular primer application may be readily determined in the manner described in H. Erlich, PCR Technology, Principles and Application for DNA Amplification (1989).
- [0024] As used herein, a "carrier" is an individual who carries a mutated or altered allele of a gene but is not affected by the disorder or disease associated with the mutation. Carriers can pass the mutation to a child or offspring in future generations, who may be affected with the disease or disorder. With respect to Fragile X syndrome, both males and females may be carriers. As used herein, the term carrier encompasses individuals having a premutation allele or a full mutation allele. Such individual carriers may also be referred to herein as "premutation carrier" (i.e., having a premutation FMR1 allele) or a "full mutation carrier" (i.e., having a full mutation FMR1 allele).
- [0025] As used herein, "5'-untranslated region (5'-UTR)" of the FMR1 gene means the region of FMR1 gene upstream of the first codon that is not translated into a protein. In one embodiment, the 5'-UTR of FMR1 gene includes promoters, enhancers and CGG repeats.
- [0026] As used herein "CGG repeat tract" means the sequence of tandem repeats of CGG trinucleotides in the 5'-UTR of the FMR1 gene. The number of CGG repeats in the CGG repeat tract may vary within an individual or between two individuals. The number of CGG repeats in the CGG repeat tract of the FMR1 gene can be at least 5, 10, 15, 20, 30, 40, 50, 75, 100, 200, 300, 500, 1,000, 2,000 or more. In some embodiments, the number of CGG repeats in the CGG repeat tract can vary from about 5 to about 2,000, to about 2,100, to about 2,300,

to about 2,500 repeats. An exemplary CGG repeat tract sequence of FMR1 gene corresponds to nucleotide positions 13833 to 13962 of SEQ ID NO: 1.

0027] As used herein "CGG repeat expansion" means a number of CGG repeats at the 5'-UTR of FMR1 gene which is greater than the normal range of 5-44 repeats. Typically normal repeat lengths are between 5-44. In some embodiments, the repeat lengths can increase to 45-54 repeats. In some embodiments, repeat lengths can increase to between 55-200 CGG repeats. Individuals having between 55-200 CGG repeats in the 5'-UTR of FMR1 gene are considered to have premutation allele. In some embodiments, the repeat lengths can increase to over 200 repeats. Individuals having greater than 200 repeats in the 5'-UTR of the FMR1 gene are considered to have a full mutation allele. The following abbreviations are used. NF: normal female; GF, intermediate female; PF, premutation female; FF; affected female; NM, normal male; GM, intermediate male; PM, premutation male; FM, affected male. "Normal" is used herein has 5-44 repeats, "Intermediate has 45-54 repeats, "Premutation" has 55-200 repeats and "Full Mutation" has >200 repeats.

[0028] As used herein, "5'-untranslated region (5'-UTR)" of the FMR2 gene means the region of FMR2 gene upstream of the first codon that is not translated into a protein. In one embodiment, the 5'-UTR of FMR2 gene includes promoters, enhancers and CCG repeats.

[0029] As used herein "CCG repeat tract" means the sequence of tandem repeats of CCG trinucleotides in the 5'-UTR of the FMR2 gene. The number of CCG repeats in the CCG repeat tract may vary within an individual or between two individuals. The number of CCG repeats in the CCG repeat tract of the FMR2 gene can be at least 5, 10, 15, 20, 30, 40, 50, 75, 100, 200, 300, 500, 1,000, 2,000 or more. In some embodiments, the number of CCG repeats in the CCG repeat tract can vary from about 5 to about 2,000, to about 2,100, to about 2,300, to about 2,500 repeats. An exemplary CCG repeat tract sequence of FMR2 gene corresponds to nucleotide positions 21 to 65 of SEQ ID NO: 75.

0030] As used herein "CCG repeat expansion" means a number of CCG repeats at the 5'-UTR of FMR2 gene which is greater than the normal range of 6-30 repeats. Typically normal repeat lengths are between 6-30. In some embodiments, the repeat lengths can increase to 31-60 repeats. In some embodiments, repeat lengths can increase to between 61-200 CCG repeats. Individuals having between 31-60 CCG repeats in the 5'-UTR of FMR2 gene are considered to have intermediate allele. Individuals having between 61-200 CCG

repeats in the 5'-UTR of FMR2 gene are considered to have premutation allele. In some embodiments, the repeat lengths can increase to over 200 repeats. Individuals having greater than 200 CCG repeats in the 5'-UTR of the FMR2 gene are considered to have a full mutation allele. The following abbreviations are used. NF: normal female; GF, intermediate female; PF, premutation female; FF; affected female; NM, normal male; GM, intermediate male; PM, premutation male; FM, affected male. "Normal" as used herein has 6-30 repeats, "Intermediate has 31-60 repeats, "Premutation" has 61-200 repeats and "Full Mutation" has >200 repeats.

[0031] As used herein, "penetrance in patients" in the context of Fragile X-associated tremor and ataxia syndrome (FXTAS) and Fragile X-associated premature ovarian insufficiency (FXPOI) means the chance of being affected by FXTAS or FXPOI as is calculated by the fraction of the population affected with FXTAS and FXPOI having premutation in FMR1 gene which may be ≥ 30 CGG repeats, ≥ 40 CGG repeats, ≥ 50 CGG repeats, ≥ 70 CGG repeats, or ≥ 100 CGG repeats. Alternatively, the permutation size may be determined relative to the population affected with FXTAS and/or FXPOI.

[0032] "Sample" or "patient sample" as used herein includes biological samples such as cells, cell lysates, tissues and body fluids. Sample may be of human or non-human origin. Sample may include, but are not limited to, amniotic fluid, biopsies, blood, blood spots, blood cells, bone marrow, amniotic fluid, cerebrospinal fluid, chorionic villi, fecal samples, excrements, fine needle biopsy samples, peritoneal fluid, plasma, pleural fluid, bronchial alveolar lavage, bronchial wash, saliva, semen, serum, sputum, tears, buccal swab, tissue, tissue homogenates, frozen tissue, paraffin sections of tissue, tissue culture media, cells, fetal cells, cell lysates, cell from culture, cell culture supernatant, fetus, embryo, and urine. A sample may include a body fluid that is "acellular."

[0033] An "acellular body fluid" is a sample that contains less than about 1% (w/w) of whole cellular material. Acellular body fluid includes samples that are naturally acellular when obtained from the individual, or may be cellular fluids that are made to be acellular. For example, plasma is an acellular body fluid derived from blood by removing whole cellular material from blood by methods known in the art (e.g., centrifugation, filtration, and the like). Serum is an acellular body fluid obtained from clotted blood after the clotted fraction is removed. Naturally occurring acellular body fluids include, for example, certain

samples of spinal fluid and urine, such as are typically obtained from a normal healthy individual.

[0034] As used herein, the term "oligonucleotide" refers to a short polymer composed of deoxyribonucleotides, ribonucleotides or any combination thereof. Oligonucleotides of the invention are generally between about 10 and about 100 nucleotides in length.

Oligonucleotides are preferably 15 to 70 nucleotides long, with 20 to 26 nucleotides being the most common. The single letter code for nucleotides is as described in the U.S. Patent Office Manual of Patent Examining Procedure, section 2422, table 1. In this regard, the nucleotide designation "R" means guanine or adenine, "Y" means thymine (uracil if RNA) or cytosine; and "M" means adenine or cytosine. An oligonucleotide may be used as a primer or as a probe.

[0035] As used herein, the term "nucleic acid" or "nucleic acid sequence" refers to an oligonucleotide, nucleotide or polynucleotide, and fragments or portions thereof, which may be single or double stranded, or partially double stranded and represent the sense or antisense strand. A nucleic acid may include DNA or RNA, and may be of natural or synthetic origin and may contain deoxyribonucleotides, ribonucleotides, or nucleotide analogs in any combination. Nucleic acid may comprise a detectable label. Although a sequence of the nucleic acids may be shown in the form of DNA, a person of ordinary skill in the art recognizes that the corresponding RNA sequence will have a similar sequence with the thymine being replaced by uracil i.e. "t" with "u".

[0036] Non-limiting examples of nucleic acid include a gene or gene fragment, genomic DNA, exons, introns, mRNA, tRNA, rRNA, ribozymes, cDNA, recombinant nucleic acid, branched nucleic acid, plasmids, vectors, isolated DNA of any sequence, isolated RNA of any sequence, synthetic nucleic acid, nucleic acid probes and primers. Nucleic acid may comprise modified nucleotides, such as methylated nucleotides and nucleotide analogs, uracyl, other sugars and linking groups such as fluororibose and thiolate, and nucleotide branches. A nucleic acid may be modified such as by conjugation, with a labeling component. Other types of modifications included in this definition are caps, substitution of one or more of the naturally occurring nucleotides with an analog, and introduction of chemical entities for attaching the polynucleotide to other molecules such as proteins, metal ions, labeling components, other nucleic acid or a solid support. Nucleic acid may include nucleic acid that has been amplified (e.g., using polymerase chain reaction).

[0037] "Genomic nucleic acid" as used herein, refers to the nucleic acid in a cell that is present in the cell chromosome(s) of an organism which contains the genes that encode the various proteins of the cells of that organism. A preferred type of genomic nucleic acid is that present in the nucleus of a eukaryotic cell. In a preferred embodiment, a genomic nucleic acid is DNA. Genomic nucleic acid can be double stranded or single stranded, or partially double stranded, or partially single stranded or a hairpin molecule. Genomic nucleic acid may be intact or fragmented (e.g., digested with restriction endonucleases or by sonication or by applying shearing force by methods known in the art). In some cases, genomic nucleic acid may include sequence from all or a portion of a single gene or from multiple genes, sequence from one or more chromosomes, or sequence from all chromosomes of a cell. As is well known, genomic nucleic acid includes gene coding regions, introns, 5' and 3' untranslated regions, 5' and 3' flanking DNA and structural segments such as telomeric and centromeric DNA, replication origins, and intergenic DNA. Genomic nucleic acid representing the total nucleic acid of the genome is referred to as "total genomic nucleic acid."

[0038] Genomic nucleic acid may be obtained by methods of extraction/purification from cellular or acellular material as is well known in the art. The ultimate source of genomic nucleic acid can be normal cells or may be cells that contain one or more mutations in the genomic nucleic acid, e.g., duplication, deletion, translocation, and transversion. Included in the meaning of genomic nucleic acid is genomic nucleic acid that has undergone recombination and may comprise rearranged genes. Also included in the meaning of genomic nucleic acid is genomic nucleic acid that has been subjected to an amplification step that increases the amount of the target sequence of interest sought to be detected relative to other nucleic acid sequences in the genomic nucleic acid. In some embodiments, genomic nucleic acid may be substantially purified from its natural environment. In another embodiment, genomic nucleic acid may not be purified from its natural environment. In some embodiments, genomic nucleic acid may be present in an environment comprising RNA, proteins such as in cell lysates, tissue homogenates.

[0039] As used herein, the term "substantially purified" in reference to oligonucleotides does not require absolute purity. Instead, it represents an indication that the sequence is relatively more pure than in the natural environment. Such oligonucleotides may be obtained by a number of methods including, for example, laboratory synthesis, restriction enzyme

digestion or PCR. A "substantially purified" oligonucleotide is preferably greater than 50% pure, more preferably at least 75% pure, and most preferably at least 95% pure.

[0040] As used herein, an oligonucleotide is "specific" for a nucleic acid if the oligonucleotide has at least 50% sequence identity with a portion of the nucleic acid when the oligonucleotide and the nucleic acid are aligned. An oligonucleotide that is specific for a nucleic acid is one that, under the appropriate hybridization and washing conditions, is capable of hybridizing to the target of interest and not substantially hybridizing to nucleic acids which are not of interest. Higher levels of sequence identity are preferred and include at least 75%, at least 80%, at least 85%, at least 90%, at least 95% and more preferably at least 98% sequence identity.

[0041] As used herein, the term "hybridize" or "specifically hybridize" refers to a process where two complementary nucleic acid strands anneal to each other under appropriately stringent conditions. Hybridizations are typically and preferably conducted with probelength nucleic acid molecules, preferably 20-100 nucleotides in length. Nucleic acid hybridization techniques are well known in the art. See, e.g., Sambrook, et al., 1989, Molecular Cloning: A Laboratory Manual, Second Edition, Cold Spring Harbor Press, Plainview, N.Y. Those skilled in the art understand how to estimate and adjust the stringency of hybridization conditions such that sequences having at least a desired level of complementarity will stably hybridize, while those having lower complementarity will not. For examples of hybridization conditions and parameters, see, e.g., Sambrook, et al., 1989, Molecular Cloning: A Laboratory Manual, Second Edition, Cold Spring Harbor Press, Plainview, N.Y.; Ausubel, F. M. et al. 1994, Current Protocols in Molecular Biology. John Wiley & Sons, Secaucus, N.J.

[0042] The term "substantially complementary" as used herein means that two sequences hybridize under stringent hybridization conditions. The skilled artisan will understand that substantially complementary sequences need not hybridize along their entire length. In particular, substantially complementary sequences comprise a contiguous sequence of bases that do not hybridize to a target or marker sequence, positioned 3' or 5' to a contiguous sequence of bases that hybridize under stringent hybridization conditions to a target or marker sequence.

[0043] The term "complement" as used herein means the complementary sequence to a nucleic acid according to standard Watson/Crick pairing rules. A complement sequence can also be a sequence of RNA complementary to the DNA sequence or its complement sequence, and can also be a cDNA.

[0044] As used herein the term "at least a portion of" in the context of nucleic acid means at least 5, at least about 15, 20, 25, 30, 35, 40, 45, 50, 75, 100, 200, 300, 400, 500, or 1,000 contiguous nucleotides or more in length.

[0045] As used herein, the term "about," unless indicated otherwise, means in quantitative terms, plus or minus 10%.

BRIEF DESCRIPTION OF THE FIGURES

[0046] FIG. 1 shows the genomic nucleic acid sequence of a portion of the human FMR1 gene including the 5'-untranslated region (SEQ ID NO: 1).

[0047] Fig. 2 A-K shows a schematic representation of the relative position of the primers for amplifying CGG repeat tract of 5'-UTR of FMR1 gene. Fig. 2A shows a forward primer (F) and reverse primer (R1) hybridizing to 5'-UTR of FMR1 gene nucleotide sequences upstream and downstream of CGG repeat tract, respectively. The reverse primers R2 and R4 , shown in Fig. 2B and 2D, respectively, hybridize to a portion of CGG repeat, the junction of 3'-end of CGG repeat tract and down stream sequences. The reverse primers R3 and R5, shown in Fig. 2C and 2E, respectively, hybridize randomly to the CGG repeat tract. Reverse primers R4 and R5, shown in Fig. 2D and 2E, respectively, further include a sequence with no homology to FMR1 gene. This sequence is complementary to primer R6 which can serve as a primer for subsequent amplification. The forward primers F1 and F3 shown in Fig. 2G and 2I, respectively, hybridize to a portion (technically the complement) of the CGG repeat and the junction of 5'-end of CGG tract and upstream sequences. The forward primers F2 and F4 shown in Fig. 2H and 2J, respectively, hybridize randomly to the CGG repeat sequence (technically the complement of this sequence). Forward primers F3 and F4 shown in Fig. 2K further include a sequence with no homology to FMR1 gene. This sequence is identical to primer F5 which can serve as a primer for subsequent amplification. After initial rounds of PCR, extended forward or reverse primers can themselves serve as both templates and primers. This results in various PCR product sizes giving the "stutter" pattern upon capillary electrophoresis.

[0048] Fig. 3 is an electropherogram showing the effect of 7-deaza-2-deoxyGTP on full mutation detection in FMR1 gene. The primers used to amplify the CGG repeat region were FMR1F, FMR1R, FMR1CCGR, and M13R. PCR master mixes were formulated with shown concentrations of 7-deaza-dGTP and dGTP. An affected female DNA sample (19/450) was used in the PCR reactions. The positions of the normal and intermediate alleles (< 55 CGG repeats) and premutation and full mutation alleles (≥ 55 CGG repeats) are indicated below the electropherograms. The arrow shows the signal for the full mutation. Replacing dGTP with 7-deaza-dGTP influences the migration and detectability of the amplified full mutation fragment.

[0049] Fig. 4 is an electropherogram showing exemplary amplification results of th CGG repeat tract of 5'-UTR of FMR1 gene obtained from normal, carrier and affected female and male fragile X patients using primers FMR1F and FMR1R and M13R in presence of 7-deaza-dGTP. The status of alleles as to number of repeats is shown for each patient. The assay distinguishes between normal females or males, including normal homozygous females, from full mutation carrier females or males. The last panel shows a sample from a mosaic male with 105/350 alleles. The 105 repeat allele is marked by an asterisk.

[0050] Fig. 5 is an electropherogram showing exemplary amplification results of the CGG repeat tract of 5'-UTR of FMR1 gene in mosaic males and females using primers FMR1F and FMR1R and M13R in presence of 7-deaza-dGTP. Samples 15, 17 from male patients, and samples 2 and 7 from female patients (identified in Table 4) have mosaic premutation and full mutations as shown. Sample 17 shows weak signal for full mutation. Other samples show robust stutter with pre- and full mutation amplification. Premutations are marked by an asterisk while mutations are marked by an arrow.

[0051] Fig. 6 is an electropherogram showing the results of assay sensitivity to detect low levels of full mutations in the 5'-UTR of FMR1 gene. DNA samples from normal males or females were mixed with affected males or females to give final concentrations of full mutation content as shown in each panel and were amplified using primers FMR1F and FMR1R and M13R in presence of 7-deaza-dGTP. The normal alleles show full scale low number repeat signal while the mutant alleles show stutter amplification. A) Set A, from mixed normal male (28 repeats) and affected male (850 repeats), shows detection of full mutation at 5% "mosaic" content. B) Set B, from mixed normal homozygous female (28/28 repeats) and affected female (19/450 repeats), shows detection of full mutation in female

DNA at 1% indicated by the arrow in the electropherogram at the right. In this instance, the full mutation allele represents only 1 of the four alleles present in the mix, and is present at very low amounts in the low "mosaic" content. C) Set C, from mixed normal male (28 repeats) and affected male (853 repeats) shows "mosaic" detection from males at 10% "mosaic" content indicated by the arrow in the electropherogram at the right. The disappearance of the normal 28 repeat allele (right most arrow in <55 repeat segment) occurs as the % full mutation is increased.

[0052] Fig. 7 is an electropherogram showing the results for detection of normal, premutation and full mutation carriers in the 5'-UTR of FMR1 gene from blood spots. Blood samples were spotted on filter cards, and punched using BSD1000 GenePunch. DNA from two spots was extracted and PCR was performed using primers FMR1F and FMR1R and M13R in presence of 7-deaza-dGTP. The number in each panel corresponds to the patient numbers in Table 5. The mutation status and CCG repeat allele content for each sample is also shown. The method can identify full mutations and distinguish normal females, including homozygous females, from carrier females.

[0053] Fig. 8 is an electropherogram showing the results of amplification of the 5'-UTR of FMR1 gene using reverse primer FMR1R (SEQ ID NO: 3), forward primer FMR1F, and M13R in presence of 7-deaza-dGTP.

[0054] Fig. 9 is an electropherogram showing the results of amplification of the 5'-UTR of FMR1 gene using reverse primer CGG1 (SEQ ID NO: 5), forward primer FMR1F, and M13R in presence of 7-deaza-dGTP. The samples used in the assay are the same as in Figure 8. The mutation status and CCG repeat allele content are shown.

[0055] Fig. 10 is an electropherogram showing the results of amplification of the 5'-UTR of FMR1 gene using reverse primer CGG2 (SEQ ID NO: 6), forward primer FMR1F, and M13R in presence of 7-deaza-dGTP. The samples used in the assay are the same as in Figure 8. The mutation status and CCG repeat allele content are shown.

[0056] Fig. 11 is an electropherogram showing the results of amplification of the 5'-UTR of FMR1 gene using reverse primer CGG3 (SEQ ID NO: 7), forward primer FMR1F, and M13R in presence of 7-deaza-dGTP. The samples used in the assay are the same as in Figure 8. The mutation status and CCG repeat allele content are shown.

[0057] Fig. 12 is an electropherogram showing the results of amplification of the 5'-UTR of FMR1 gene using reverse primer CGG4 (SEQ ID NO: 8), forward primer FMR1F, and M13R in presence of 7-deaza-dGTP. The samples used in the assay are the same as in Figure 8. The mutation status and CCG repeat allele content are shown.

[0058] Fig. 13 is an electropherogram showing the results of amplification of the 5'-UTR of FMR1 gene using reverse primer CGG3-1 (SEQ ID NO: 9), forward primer FMR1F, and M13R in presence of 7-deaza-dGTP. The samples used in the assay are the same as in Figure 8.

[0059] Fig. 14 is an electropherogram showing the results of amplification of the 5'-UTR of FMR1 gene using reverse primer CGG 4-1 (SEQ ID NO: 10), forward primer FMR1F, and M13R in presence of 7-deaza-dGTP. The samples used in the assay are the same as in Figure 8.

[0060] Fig. 15 shows exemplary arrangement of CCG triplet sequences relative to each other in the primer.

[0061] Fig. 16 shows exemplary arrangement of CGG triplet sequences relative to each other in the primer.

[0062] Fig. 17 shows an exemplary testing algorithm for samples submitted for Fragile X testing. Samples received will be tested using two different PCR based assays. Reaction 1 can be used for sizing normal, intermediate and small premutations, as well as gender determination using Capillary Southern Assay (Strom et al. Genet Med. 2007; 9(4):199-207). Reaction 2 (triplet primed amplification of the invention) will detect the presence of expanded alleles and full mutations. Samples showing evidence of stutter or full mutation amplification will be further analyzed by Southern blot analysis for methylation status and sizing. As shown in the figure, the use of triplet primed amplification will significantly reduce the number of Southern blot analyses to be performed because the method clearly distinguishes normal homozygous females, which represent ~ 95% of all female samples from full mutation carriers.

[0063] Fig. 18 shows an exemplary result of an CCG repeat-primed PCR assay and subsequent confirmation by Capillary Southern Analysis (CSA) for the detection of CGG repeats in the 5'-UTR of FMR1 gene. Figure 18A is an electropherogram showing the results

CCG repeat-primed PCR assay using reverse primer FMR1R (SEQ ID NO: 3), forward primer FMR1F, and M13R in presence of 7-deaza-dGTP. Fig 18B shows the reconfirmation of the results CCG repeat-primed PCR assay using the CSA method with gender confirmation.

[0064] Fig. 19 shows the distribution of the 48 premutation alleles of the FMR1 gene identified in DNA samples submitted for Cystic Fibrosis screening.

[0065] Fig. 20 shows the identification of a full mutation allele (> 200 CGG repeats) of the FMR1 gene in a sample obtained from a Caucasian female. Fig 20A is an electropherogram showing a late migrating fragment on capillary electrophoresis as an evidence of full mutation in the FMR1 gene. Fig. 20B shows the size separation of restriction endonuclease digested fragments to indicate the methylation status of CGG repeat tract of the 5'-UTR of FMR1 gene. The appearance of additional DNA fragments (shown by arrow) along with the expected fragment sizes is indicative of heterogeneity in DNA methylation.

[0066] Fig. 21 shows an exemplary nucleic acid sequence of the 5'-UTR of FMR2 gene comprising the CCG repeats (SEQ ID NO: 75).

Fig. 22 A-K shows a schematic representation of the relative position of the amplification primers for amplifying CCG repeat tract in the 5'-UTR of FMR2 gene. Fig. 22A shows a forward primer (F) and reverse primer (R1) hybridizing to 5'-UTR of FMR2 gene nucleotide sequences upstream and downstream of CCG repeat tract respectively. The reverse primers R2 and R4, shown in Fig. 22B and 22D, respectively, hybridize to a portion of CCG repeat, the junction of 3'-end of CCG repeat tract and down stream sequences. The reverse primers R3 and R5, shown in Fig. 22C and 22E, respectively, hybridize randomly to the CCG repeat tract. Reverse primers R4 and R5, shown in Fig. 22D and 22E, respectively, further include a sequence with no homology to FMR2 gene. This sequence is complementary to primer R6 which can serve as a primer for subsequent amplification. The forward primers F1 and F3 shown in Fig. 22G and 22I, respectively, hybridize to a portion (technically the complement) of the CCG repeat and the junction of 5'-end of CCG tract and upstream sequences. The forward primers F2 and F4 shown in Fig. 22H and 22J, respectively, hybridize randomly to the CCG repeat sequence (technically the complement of this sequence). Forward primers F3 and F4 shown in Fig. 22K further include a sequence with no homology to FMR2 gene. This sequence is identical to primer F5 which can serve as

a primer for subsequent amplification. After initial rounds of PCR, extended forward or reverse primers can themselves serve as both templates and primers. This results in various PCR product sizes giving the "stutter" pattern upon capillary electrophoresis.

DETAILED DESCRIPTION OF THE INVENTION

[0068] Repeat expansion disease includes any of about two dozen human diseases displaying Mendelian inheritance patterns shown to be caused by expansions of intrinsically polymorphic tandem repeats, mainly involving different trinucleotide motifs but also longer repetitive sequences up to 12-mers (Table 1). A characteristic of an allele containing an expanded tandem repeat is an excessive instability in successive generations (dynamic mutations). Furthermore, these alleles can differ in lengths among cell populations of the same organism (mosaicism). One type of repeat expansion disease is the trinucleotide repeat disorders (e.g., fragile X syndrome, myotonic dystrophy 1, etc.), the most abundant form of repeat expansion diseases. These diseases exhibit intergenerational repeat instability with a tendency towards further expansion of the tandem repeat. Increased repeat lengths in successive generations can lead to an earlier age of onset in affected individuals and/or an accentuation of clinical symptoms. The methods of the present invention which are experimentally exemplified for determining an expansion of the CGG repeat tract of the FMR1 gene also are applicable to detecting the same repeat tracts in genes associated with other diseases as shown in Table 1.

Table 1. Exemplary tandem repeat expansion diseases

Disease (gene designation)	Protein	Tandem repeat sequence	Repeat Localization	Normal repeat range	Expanded repeat
Dentatorubral- pallidoluysian atrophy (DRPLA)	atrophin-1	CAG	coding	6-35	49-88
Progressive myclonus epilepsy (EPM1)#	CSTB	CCCCGC CCCGCG	5' UTR	2-3	50-80
Fragile XA- syndrome (FRAXA)*	FMR1	CGG	5' UTR	5-44	>230
Fragile XE- syndrome (FMR1)*	FMR-2	CCG	5' UTR	6-30	>200

Friedreich ataxia (FRDA)#	Frataxin	GAA	intronic	7-34	>100
Huntington disease (huntingtin)	huntingtin	CAG	coding	4-35	36-250
Huntington disease-like 2 (JPH-3)	junctophilin- 3	CTG	coding	6-27	>40-60
Spinobulbar musc. atrophy, Kennedy disease (AR)*	androgen rec.	CAG	coding	9-36	38-62
Myotonic dystrophy 1 (DM1)	DMPK/SIX5	CTG	3' UTR	5-37	>50
Myotonic dystrophy 2 (DM2)	ZNF9	ССТG	intronic	10-26	75->11000
Spinocerebellar ataxia 1 (SCA1)	ataxin-1	CAG	coding	6-44	39-82
Spinocerebellar ataxia 2 (SCA2)	ataxin-2	CAG	coding	15-31	36-63
Spinocerebellar ataxia 3, Machado-Joseph disease (SCA3)	Ataxin-3	CAG	coding	12-40	55-84
Spinocerebellar ataxia 6 (SCA6)	CACNA1A	CAG	coding	4-18	21-33
Spinocerebellar ataxia 7 (SCA7)	Ataxin-7	CAG	coding	4-35	37-306
Spinocerebellar ataxia 8 (SCA8)	not known	СТG	3' UTR	16-37	110-250
Spinocerebellar ataxia 10 (SCA10)	Ataxin-10	ATTCT	intronic	10-22	800-4600
Spinocerebellar ataxia 12 (SCA12)	PP2A-PR55B	CAG	promotor	7-28	66-78
Spinocerebellar ataxia 17 (SCA17)	ТВР	CAG	coding	25-42	47-63
Oculopharyngeal muscular dystrophy (PABPN1)	PABPN1	GCG	coding	6	7-13
Synpolydactyly, type II (HOXD13)	HOXD13	GCN	coding	15	22-24

^{* =} X chromosome; # = autosomal recessive; undesignated = autosomal dominant

[0069] Methods of detecting a particular nucleic acid segment of interest in a sample of nucleic acids are provided. In particular embodiments, the particular nucleic acid segment of interest is a tandem repeat and the method is used to determine information about the size of such tandem repeat. In some embodiments, the methods are useful for detecting CGG repeats in the FMR1 gene. In some embodiments, the methods are useful for detecting CCG repeats in the FMR2 gene. The methods are useful in detecting over 2,000 CGG and/or CCG repeats in the highly GC rich, CGG containing 5'-UTR the region of FMR1 gene and CCG containing 5'-UTR the region of FMR1 gene respectively. The methods involve amplifying the CGG containing 5'-UTR region of the FMR1 gene and/or CCG containing 5'-UTR region of the FMR2 gene followed by detecting the amplicons. The methods allow detection of CGG expansions and full mutations in FMR1 gene of both males and females. The methods also allow detection of CCG expansions and full mutations in FMR2 gene of both males and females. The methods can be used on small sample amounts such as a blood spot that are typically available in newborn screening programs. The methods are useful to determine the frequency of Fragile X syndrome premutations and full mutations in nonselected unbiased populations undergoing routine carrier screening.

[0070] FMR1 Gene

[0071] Human FMR1 gene is located in chromosome X. Fragile X syndrome (FXS) is caused by expansion and subsequent methylation of the CGG trinucleotide repeat in the *FMR1* 5'-untranslated region (5'-UTR). Typically, normal repeat lengths are between 5-44, intermediate are between 45-54, premutation are between 55-200 CGG, while full mutation are greater than 200. Exemplary FMR1 genomic nucleic acid sequence comprising the 5'-UTR include, but is not limited to, sequence in GenBank accession number NG_007529 and in L29074. These sequences are incorporated herein by reference. Exemplary FMR1 nucleic acid sequence comprising the 5'-UTR, as well as the start codon and a portion of the beginning of the coding sequence as set forth in GenBank accession number L29074 is disclosed as SEQ ID NO: 12 in US Patent Application Publication 2008/0113355. The sequence of SEQ ID NO: 12 of US 2008/0113355 is incorporated herein by reference. The CGG repeat tract starts at position 13833 and ends in position 13962 and is depicted herein in SEQ ID NO: 1.

[0072] FMR2 Gene

[0073] Human FMR2 gene is located in chromosome Xq28, 600 kb distal to FMR1 gene. The FMR2 (FRAXE) mutation is an expansion of the CCG containing repeats in the 5'-UTR of FMR2 gene that result in the methylation of a nearby CpG island. Typically, normal CCG repeat lengths are between 6-30, intermediate are between 31-60, permutations are between 61-200 and full mutations are greater than 200. Exemplary FMR1 genomic nucleic acid sequence comprising the 5'-UTR include, but is not limited to, sequence in GenBank accession number AH005569. This sequence is incorporated herein by reference. An exemplary sequence of the 5'-UTR of FMR2 gene comprising the CCG repeats is shown in Fig. 21 and depicted as SEQ ID NO: 75.

[0074] Biological Sample Collection and Preparation

[0075] The biological sample may be obtained from a stage of life such as a fetus, young adult, adult, and the like. In one example, subjects are humans being tested for the existence expanded CGG repeats in the 5'-UTR of FMR1 gene. In another example, subjects are humans being tested for the existence expanded CCG repeats in the 5'-UTR of FMR2 gene. In one embodiment, the sample may be obtained from an individual who is suspected of having a fragile X syndrome. In another embodiment, the sample may be obtained from a healthy individual who is assumed of having no disease, or a genetic abnormality. In another embodiment, the sample may be obtained from an individual suspected of being a carrier of fragile X syndrome but not affected by the disease.

[0076] <u>Sample Collection:</u> Methods of obtaining samples are well known to those of skill in the art and include, but are not limited to, aspirations, tissue sections, drawing of blood or other fluids, surgical or needle biopsies, collection of paraffin embedded tissue, collection of body fluids, collection of stool, and the like.

[0077] The methods can be used to assist with prenatal diagnosis of fragile X or other triplet repeat disease using any type of embryonic or fetal cell or nucleic acid containing body fluid. Fetal cells can be obtained through the pregnant female, or from a sample of an embryo. Thus, fetal cells are present in amniotic fluid obtained by amniocentesis, chorionic villi aspirated by syringe, percutaneous umbilical blood, a fetal skin biopsy, a blastomere from a four-cell to eight-cell stage embryo (pre-implantation), or a trophectoderm sample from a blastocyst (pre-implantation or by uterine lavage).

[0078] Genomic nucleic acid: In one embodiment, nucleic acid is genomic nucleic acid. Genomic nucleic acid may be intact. In another embodiment, genomic nucleic acid may be fragmented (e.g., digested with restriction endonucleases, or by sonication or by applying shearing force by methods known in the art).

[0079] Sample Preparation: The nucleic acid (DNA or RNA) may be isolated from the sample according to any methods well known to those of skill in the art. If necessary, the sample may be collected or concentrated by centrifugation and the like. The sample may be subjected to lysis, such as by treatments with enzymes, heat, surfactants, ultrasonication, or a combination thereof. The lysis treatment is performed in order to obtain a sufficient amount of nucleic acid. The sample may be subjected to liquid chromatography to partially purify the genomic nucleic acid.

[0080] Suitable DNA isolation methods include phenol and chloroform extraction. See Maniatis et al., Molecular Cloning, A Laboratory Manual, 2d, Cold Spring Harbor Laboratory Press, page 16.54 (1989). Numerous commercial kits also yield suitable DNA including, but not limited to, QIAampTM mini blood kit, Agencourt GenfindTM, Roche Cobas® Roche MagNA Pure® or phenol:chloroform extraction using Eppendorf Phase Lock Gels®.

[0081] Total DNA (e.g., genomic, mitochondrial, microbial, viral,) can be purified from any biological sample using commercially available kits e.g. QIAamp DNA and QIAamp DNA Blood mini kits, Qiagen M96 robot and reagents, Qiagen Gentra robot and reagents, and Qiagen 9604 reagents (Qiagen, Valencia, CA). In one example, blood can be spotted on Guthrie cards. Blood spots can be punched from each card using BSD1000 GenePunch Instrument and DNA was extracted using Qiagen BioSprint reagents.

[0082] Genomic DNA may be isolated from cells or tissues using standard methods, see, *e.g.*, Sambrook, et al., 1989, Molecular Cloning: A Laboratory Manual, Second Edition, Cold Spring Harbor Press, Plainview, NY.

[0083] In another embodiment, sample nucleic acid may be mRNA or cDNA generated from mRNA or total RNA. RNA can be isolated from cells or tissue samples using standard techniques, see, *e.g.*, Sambrook, et al., 1989, Molecular Cloning: A Laboratory Manual, Second Edition, Cold Spring Harbor Press, Plainview, NY. In addition, kits for isolating mRNA and synthesizing cDNA are commercially available e.g. RNeasy Protect Mini kit, RNeasy Protect Cell Mini kit from Qiagen.

[0084] Nucleic acid need not be extracted, but may be made available by suitable treatment of cells or tissue such as described in US Patent 7,521,213.

[0085] Nucleic Acid Amplification

[0086] Nucleic acid comprising triplet repeat sequences can be amplified by various methods known to the skilled artisan. In one embodiment, nucleic acid comprising CGG repeat sequences of the FMR1 gene or nucleic acid comprising CCG repeat sequences of the FMR 2 gene can be amplified by PCR. PCR methods to amplify CGG repeats and CCG repeats include, but are not limited to, nested PCR which involves a set of outer primer and a set of inner primer, semi-nested PCR which involves a common primer, an outward primer and an inner primer, and standard PCR with a primer pair. Efficient amplification of the GC rich CGG repeat and CCG repeat containing 5-UTR region of FMR 1 and FMR 2 respectively will depend on several parameters including, but not limited to, primer location, primer sequence, nucleotide concentrations, use of a dGTP analogue such as 7-deaza-2-deoxyGTP, and inclusion of additives to increase the yield, specificity and consistency.

[0087] Primer locations and primer sequences:

[0088] Primer locations for the detection of CGG repeats in FMR1 gene

[0089] In one embodiment, the forward and reverse primers may hybridize to regions in the 5'-UTR of the FMR1 gene outside the CGG repeat sequence as shown in Figure 2A. In another embodiment, the reverse primer may comprise one or more CCG triplet sequences and include sequence complementary to the junction sequence of the 3'-end of the CGG repeat tract and flanking sequence and a portion of the sequence downstream of the CGG repeat tract as shown in Figure 2B. One or more CCG triplets will hybridize to the CGG repeats and permit amplification of the CGG repeat sequences. In yet another embodiment, the reverse primer may comprise one or more CCG triplet sequences without any flanking and junction sequences such that the reverse primers may hybridize randomly to the CGG repeat sequence and initiate DNA synthesis to generate a broad spectrum of amplicons of varying length or stutter. Exemplary reverse primers are shown in Fig. 2A-F.

[0090] In another embodiment, the forward primer may comprise one or more CGG triplet sequences, the junction sequence of the 5'-end of the CGG repeat tract and flanking sequence and a portion of the sequence upstream of the CGG repeat tract as shown in Figure 2G. One or more CGG triplets will hybridize to the CCG repeats of the complementary strand and

permit amplification of the CGG repeat sequences. In yet another embodiment, the forward primer may comprise one or more CGG triplet sequences without any flanking and junction sequences such that the forward primers may hybridize randomly to the complementary strand of CGG repeat sequence and initiate DNA synthesis to generate a broad spectrum of amplicons of varying length or stutter. Exemplary forward primers are shown in Fig. 2G-K.

[0091] Primer locations for the detection of CCG repeats in FMR2 gene

[0092] In one embodiment, the forward and reverse primers may hybridize to regions in the 5'-UTR of the FMR2 gene outside the CCG repeat sequence as shown in Figure 22A. In another embodiment, the reverse primer may comprise one or more CGG triplet sequences and include sequence complementary to the junction sequence of the 3'-end of the CCG repeat tract and flanking sequence and a portion of the sequence downstream of the CCG repeat tract as shown in Figure 22B. One or more CGG triplets will hybridize to the CCG repeats and permit amplification of the CCG repeat sequences. In yet another embodiment, the reverse primer may comprise one or more CGG triplet sequences without any flanking and junction sequences such that the reverse primers may hybridize randomly to the CCG repeat sequence and initiate DNA synthesis to generate a broad spectrum of amplicons of varying length or stutter. Exemplary reverse primers are shown in Fig. 22A-F.

[0093] In another embodiment, the forward primer may comprise one or more CCG triplet sequences, the junction sequence of the 5'-end of the CCG repeat tract and flanking sequence and a portion of the sequence upstream of the CCG repeat tract as shown in Figure 22G. One or more CCG triplets will hybridize to the CGG repeats of the complementary strand and permit amplification of the CCG repeat sequences. In yet another embodiment, the forward primer may comprise one or more CCG triplet sequences without any flanking and junction sequences such that the forward primers may hybridize randomly to the complementary strand of CCG repeat sequence and initiate DNA synthesis to generate a broad spectrum of amplicons of varying length or stutter. Exemplary forward primers are shown in Fig. 22G-K.

[0094] In some embodiments, the one or more primers can be labeled with a detectable label. Exemplary detectable label includes but not limited to fluorophores, radiolabels, chemiluminiscent compound, colorimetric labels (e.g., colloidal gold), magnetic labels (e.g., DynabeadsTM), biotin, digoxigenin, haptens, proteins, enzymes.

[0095] In some embodiments, the forward or the reverse or both primers may comprise degenerate bases. In some embodiments, the forward or the reverse or both primers may comprise universal bases. Exemplary universal bases include 3-nitropyrrole, 5'-nitroindole.

[0096] In some embodiments, the forward or the reverse or both primers may further comprise an optional sequence at the 5'-end that has no homology with the target sequence of SEQ ID NO: 1 or SEQ ID NO: 75 (see, Figures 2 D-F, 2 I-K, 22 D-F, and 22 I-K). Such 5'-end sequences can be 6-50 nucleotides in length, 10-30 nucleotides in length, or 15-25 nucleotides in length. In one embodiment, the optional 5'-end sequence may be from bacteriophage M13 or from T7 promoter. In another embodiment, the optional 5'-end sequence may be a random sequence. After a few rounds of amplification, the resulting amplicons will comprise the optional 5'-end sequence. These 5'-end sequences can be utilized as template for further round of amplification using primers hybridizing to these 5'-end sequences. Exemplary M13 5'-end sequence is shown below and listed as SEQ ID NO: 11.

CAGGAAACAGCTATGAC (SEQ ID NO: 11)

Exemplary T7 promoter sequence is shown below and listed as SEQ ID NO: 70.

TAATACGACTCACTATAGGG (SEQ ID NO: 70)

[0097] In some embodiments, the PCR reaction may be a semi-nested PCR, such that there is a common forward primer and an inner and outer reverse primer as shown in Fig. 2F and 22F.

Primers flanking the CGG repeat region of FMR1 gene: Methods for designing suitable primers that hybridize to the 5'-UTR of the FMR1 gene as exemplified in SEQ ID NO: 1 (or SEQ ID NO: 12 of US 2008/0113355) are known in the art such as H. Erlich, PCR Technology, Principles and Application for DNA Amplification (1989). In one embodiment, the 5'-primer (forward primer) can hybridize to the 5'-UTR within 500 or fewer nucleotides, within 300 or fewer nucleotides, within 200 or fewer nucleotides, within 100 or fewer nucleotide position of 13833 of SEQ ID NO: 1. In another embodiment, the 3'-primer (reverse primer) can hybridize to the 5'-UTR within 500 or fewer nucleotides, within 300 or fewer nucleotides, within 300 or fewer nucleotides, within 50 or

fewer nucleotides, within 25 or fewer nucleotides of 3' of nucleotide position of 13962 of SEQ ID NO: 1. Exemplary primers hybridizing to 5'-UTR flanking the CGG repeat region of FMR1 gene are provided below:

Forward: 5'-GCTCAGCTCCGTTTCGGTTTCACTTCCGGT-3' (SEQ ID NO: 12)

Reverse: 5-AGCCCCGCACTTCCACTTTCGGTTTCACTTCCGGT-3' (SEQ ID NO: 13)

[0099] Primers flanking the CCG repeat region of FMR2 gene: Methods for designing suitable primers that hybridize to the 5'-UTR of the FMR2 gene as exemplified in SEQ ID NO: 75 are known in the art such as H. Erlich, PCR Technology, Principles and Application for DNA Amplification (1989). In one embodiment, the 5'-primer (forward primer) can hybridize to the 5'-UTR within 500 or fewer nucleotides, within 300 or fewer nucleotides, within 200 or fewer nucleotides, within 100 or fewer nucleotides, within 50 or fewer nucleotides, within 25 or fewer nucleotides 5' of nucleotide position of 21 of SEQ ID NO: 75. In another embodiment, the 3'-primer (reverse primer) can hybridize to the 5'-UTR within 500 or fewer nucleotides, within 300 or fewer nucleotides, within 200 or fewer nucleotides, within 20 or fewer nucleotides, within 25 or fewer nucleotides of 3' of nucleotide position of 65 of SEQ ID NO: 75. Exemplary primers hybridizing to 5'-UTR flanking the CCG repeat region of FMR2 gene are provided below:

Forward: 5'- CGCCGCCTGTGCAGCCGCTG -3' (SEQ ID NO: 76)

Reverse: 5- GCCCCGGCTGCCGCGCGCG-3' (SEQ ID NO: 77)

[0100] Reverse primers with CCG triplet sequences for detecting CGG repeat tract in the 5'-UTR of FMR1 gene: Reverse primers with one or more CCG triplet sequences will hybridize to the CGG repeats sequences and permit amplification of the repeat region. The reverse primers may vary in the number of CCG triplet sequences from 1-40 CCG triplets. In some embodiments, the primers may comprise at least: 2, 3, 4, 5, 6, 7, 8, 9, 10, 11, 12, 13, 14, 15 or more CCG triplets. In some embodiments, the primers may comprise at least 10 consecutive CCG triplets. In some embodiments, the primers may comprise less than three CCG triplets. In some embodiment, the primers may comprise less than three consecutive CCG triplets. In one embodiment, the primers comprise two consecutive CCG triplets. In some embodiments, the primers may comprise at least two sets of two consecutive CCG triplets are

separated from each other by at least one nucleotide. In some embodiments, the linking nucleotides can be A or T. In some embodiments, the nucleotides linking the two consecutive CCG triplets (CCG)₂ may form a stem loop structure. In some embodiments, the primers may comprise two consecutive CCG triplets (CCG)₂ and additional nucleotides at the 5'- and 3'- end. In some embodiments, two sets of two consecutive CCG triplets (CCG)₂ are linked to each other through non-phosphodiester linkage. Exemplary non-nucleotide linkers include C4, C6, C8, C10 linkers, polyethylene glycol (PEG) linkers: PEG₄, PEG₆, PEG₁₀, PEG₁₂, peptide nucleic acid linkage, phosphorothoioester linkage. Exemplary arrangements of CCG triplets relative to each other are shown in Figure 15. Optionally the primers comprising one or more CCG triplets of Figure 15 may further comprise a sequence at the 5'-end that has at least 90% homology to the sequence directly downstream (i.e. at the 3'-end) of the CGG repeat tract. In some embodiments, such sequence can be at least 3, 5, 10 or more nucleotides.

[0101] In some embodiments, reverse primers comprising one or more CCG triplets may further comprise an optional sequence at the 5'-end that has no homology with the target sequence of SEQ ID NO: 1. In some embodiments, the optional sequence at the 5'-end that has no homology to human genome. In some embodiments, reverse primers comprising one or more CCG triplets may comprise a second portion and a third portion at the 5'-end of CCG triplets segment such that the second portion comprises a sequence that has at least 90% homology to the sequence directly downstream (i.e. at the 3'-end) of the CGG repeat tract and the third portion, which is at the 5'-end of the second portion, comprises a sequence that has no homology with the target sequence of SEQ ID NO: 1 (or any of the FMR1 gene).

[0102] Exemplary reverse primer sequences comprising CCG triplets are shown in the table below.

Table 2: Reverse primers having one or more CCG triplets

Reverse Primer Sequence	SEQ ID NO:
CGCCGCCCC	14
CCGCCGACCGCCG	15
CCGCCGTCCGCCG	16
CCGCCGCCCG	17
CCGCCGCCGC	18
CGCCGCCGACGCCGCC	19
CGCCGCCGCCGCC	20

Reverse Primer Sequence	SEQ ID
-	NO:
CGCCGCCGCCGCCC	21
CGCCGCCGCCGCC	22
CTCGAGGCCCAGCCGCCG	23
CTCGAGGCCCAGCCGTCGCCG	24
CTCGAGGCCCAGCCGACGCCG	25
CTCGAGGCCCAGCCGGCGCCG	26
CAGGAAACAGCTATGACCCGCCG	27
CAGGAAACAGCTATGACGCCGCCC	28
CAGGAAACAGCTATGACCCGCCGACCGCCG	29
CAGGAAACAGCTATGACCCGCCGTCCGCCG	30
CAGGAAACAGCTATGACCCGCCGGCCGCCG	31
CAGGAAACAGCTATGACCCGCCGCCGCCG	32
CAGGAAACAGCTATGACGCCGCCGACGCCGCCC	33
CAGGAAACAGCTATGACGCCGCCGTCGCCGCCC	34
CAGGAAACAGCTATGACGCCGCCGCCGCCCCCCC	35
CAGGAAACAGCTATGACGCCGCCGCCGCCGCC	36
CAGGAAACAGCTATGACCTCGAGGCCCACCGCCG	37
CAGGAAACAGCTATGACCTCGAGGCCCAGCCGTCGCCGCCG	38
CAGGAAACAGCTATGACCTCGAGGCCCAGCCGACGCCGCCG	39
CAGGAAACAGCTATGACCTCGAGGCCCAGCCGGCGCCCCG	40
CAGGAAACAGCTATGACCCTCGAGGCCCAGCCGCCGCCG	71
CC	
CAGGAAACAGCTATGACCCTCGAGGCCCAGCCGCCGCCGCCGCCGCCCGC	72
CAGGAAACAGCTATGACCCTCGAGGCCCAGCCGCCCC	73
CAGGAAACAGCTATGACCCTCGAGGCCCAGCCGCCG	74

[0103] Forward primers with CGG triplet sequences for detecting CGG repeat tract in the 5'-UTR of FMR1 gene: Forward primers with one or more CGG triplet sequences will hybridize to the complementary strand of CGG repeats sequences and permit amplification of the repeat region. The forward primers may vary in the number of CGG triplet sequences from 1-40 CGG triplets. In some embodiments, the primers may comprise at least: 2, 3, 4, 5, 6, 7, 8, 9, 10, 11, 12, 13, 14, 15 or more CGG triplets. In some embodiments, the primers may comprise at least 10 consecutive CGG triplets. In some embodiments, the primers may comprise less than three CGG triplets. In some embodiments, the primers may comprise less than three consecutive CGG triplets. In one embodiment, the primers comprise two consecutive CGG triplets. In some embodiments, the primers may comprise at least two sets of two consecutive CGG triplets: (CGG)₂. In some embodiments, at least two sets of two consecutive CGG triplets are separated from each other by at least one nucleotide. In some embodiments, the linking nucleotides can be A or T. In some embodiments, the nucleotides

linking the two consecutive CGG triplets (CGG)₂ may form a stem loop structure. In some embodiments, the primers may comprise two consecutive CGG triplets (CGG)₂ and additional nucleotides at the 5'- and 3'- end. In some embodiments, two sets of two consecutive CGG triplets (CGG)₂ are linked to each other through non-phosphodiester linkage. Exemplary non-nucleotide linkers include C4, C6, C8, C10 linkers, polyethylene glycol (PEG) linkers: PEG₄, PEG₆, PEG₁₀, PEG₁₂, peptide nucleic acid linkage, phosphorothoioester linkage. Exemplary arrangements of CGG triplets are shown in Figure 16. Optionally the primers comprising one or more CGG triplets of Figure 16 may further comprise a sequence at the 5'-end that has at least 90% homology to the sequence directly upstream (i.e. at the 5'-end) of the CGG repeat tract. In some embodiments, such sequence can be at least 3, 5, 10 or more nucleotides. In some embodiments, forward primers comprising one or more CGG triplets of Figure 16 may further comprise an optional sequence at the 5'-end that has no homology with the target sequence of SEQ ID NO: 1. In some embodiments, the optional sequence at the 5'-end that has no homology to human genome. In some embodiments, forward primers comprising one or more CGG triplets of Figure 16 may comprise a second portion and a third portion at the 5'-end of one or more CGG triplets such that the second portion comprises a sequence that has at least 90% homology to sequence directly upstream (i.e. at the 5'-end) of the CGG repeat tract, and the third portion, which is at the 5'-end of the second portion, comprises a sequence that has no homology with the target sequence of SEQ ID NO: 1 (or any of the FMR1 gene).

[0104] Exemplary forward primer sequences comprising CCG triplets are shown in the table below.

Table 3: Forward primers having one or more CGG triplets

Forward Primer Sequence	SEQ ID NO:
CGCGGCGCG	41
CGGCGGACGGCGG	42
CGGCGGTCGGCGG	43
CGGCGGCGGCGG	44
CGGCGGCCGG	45
GGCGGCGAGGCGGCG	46
GGCGGCGGCGCG	47
GGCGGCGGGGGCGCG	48
GGCGGCGGCGGCG	49
CGTGCIGCAGCGCGG	50

Forward Primer Sequence	SEQ ID NO:
CGTGCIGCAGCGCGGTGGCGG	51
CGTGCIGCAGCGCGGAGGCGG	52
CGTGCIGCAGCGGGGGGGG	53
CAGGAAACAGCTATGACCGGCGG	54
CAGGAAACAGCTATGACCGGCGCG	55
CAGGAAACAGCTATGACCGGCGGACGGCGG	56
CAGGAAACAGCTATGACCGGCGGTCGGCGG	57
CAGGAAACAGCTATGACCGGCGGCGGCGG	58
CAGGAAACAGCTATGACCGGCGGCGGCGG	59
CAGGAAACAGCTATGACCCGGCGGAGGCGGCGCG	60
CAGGAAACAGCTATGACCCGGCGGTGGCGGCGCG	61
CAGGAAACAGCTATGACCCGGCGGGGGGGGGGGGG	62
CAGGAAACAGCTATGACCCGGCGGCGGCGGCG	63
CAGGAAACAGCTATGACCCGGCGGTCGGCGGCG	64
CAGGAAACAGCTATGACCCGGCGGCGGCGGCG	65
CAGGAAACAGCTATGACCCGGCGGACGGCGGCG	66
CAGGAAACAGCTATGACCCGTGCIGCAGCGCGGTGGCGGCGG	67
CAGGAAACAGCTATGACCCGTGCIGCAGCGCGGAGGCGGCGG	68
CAGGAAACAGCTATGACCCGTGCIGCAGCGCGGGGGGGGGG	69

[0105] Reverse primers with CGG triplet sequences for detecting CCG repeat tract in the 5'-UTR of FMR2 gene: Reverse primers with one or more CGG triplet sequences will hybridize to the CCG repeats sequences and permit amplification of the repeat region. The reverse primers may vary in the number of CGG triplet sequences from 1-40 CGG triplets. In some embodiments, the primers may comprise at least: 2, 3, 4, 5, 6, 7, 8, 9, 10, 11, 12, 13, 14, 15 or more CGG triplets. In some embodiments, the primers may comprise at least 10 consecutive CGG triplets. In some embodiments, the primers may comprise less than three CGG triplets. In some embodiments, the primers may comprise less than three consecutive CGG triplets. In one embodiment, the primers comprise two consecutive CGG triplets. In some embodiments, the primers may comprise at least two sets of two consecutive CGG triplets: (CGG)₂. In some embodiments, at least two sets of two consecutive CGG triplets are separated from each other by at least one nucleotide. In some embodiments, the linking nucleotides can be A or T. In some embodiments, the nucleotides linking the two consecutive CGG triplets (CGG)₂ may form a stem loop structure. In some embodiments, the primers may comprise two consecutive CGG triplets (CGG)₂ and additional nucleotides at the 5'- and 3'- end. In some embodiments, two sets of two consecutive CGG triplets (CGG)₂ are linked to each other through non-phosphodiester linkage. Exemplary non-nucleotide linkers

include C4, C6, C8, C10 linkers, polyethylene glycol (PEG) linkers: PEG₄, PEG₆, PEG₁₀, PEG₁₂, peptide nucleic acid linkage, phosphorothoioester linkage. Exemplary arrangements of CGG triplets are shown in Figure 16. Optionally the primers comprising one or more CGG triplets of Figure 16 may further comprise a sequence at the 5'-end that has at least 90% homology to the sequence directly downstream (i.e. at the 3'-end) of the CCG repeat tract. In some embodiments, such sequence can be at least 3, 5, 10 or more nucleotides. In some embodiments, reverse primers comprising one or more CGG triplets of Figure 16 may further comprise an optional sequence at the 5'-end that has no homology with the target sequence of SEQ ID NO: 75. In some embodiments, the optional sequence at the 5'-end that has no homology to human genome. In some embodiments, reverse primers comprising one or more CGG triplets of Figure 16 may comprise a second portion and a third portion at the 5'-end of one or more CGG triplets such that the second portion comprises a sequence that has at least 90% homology to sequence directly downstream (i.e. at the 3'-end) of the CCG repeat tract, and the third portion, which is at the 5'-end of the second portion, comprises a sequence that has no homology with the target sequence of SEQ ID NO: 75 (or any of the FMR2 gene).

[0106] Exemplary reverse primer sequences comprising CGG triplets for detecting CCG repeat tract in the 5'-UTR of FMR2 gene can be any of the CGG triplet containing forward primers used to detect CGG repeat tract in the 5'-UTR of FMR1 gene as listed in Table 3 above.

[0107] Forward primers with CCG triplet sequences for detecting CCG repeat tract in the 5'-UTR of FMR2 gene: Forward primers with one or more CCG triplet sequences will hybridize to the complementary strand of CCG repeats sequences and permit amplification of the repeat region. The forward primers may vary in the number of CCG triplet sequences from 1-40 CCG triplets. In some embodiments, the primers may comprise at least: 2, 3, 4, 5, 6, 7, 8, 9, 10, 11, 12, 13, 14, 15 or more CCG triplets. In some embodiments, the primers may comprise at least 10 consecutive CCG triplets. In some embodiments, the primers may comprise less than three CCG triplets. In some embodiments, the primers may comprise less than three consecutive CCG triplets. In one embodiment, the primers comprise two consecutive CCG triplets. In some embodiments, the primers may comprise at least two sets of two consecutive CCG triplets: (CCG)₂. In some embodiments, at least two sets of two consecutive CCG triplets are separated from each other by at least one nucleotide. In some embodiments, the linking nucleotides can be A or T. In some embodiments, the nucleotides

linking the two consecutive CCG triplets (CCG)₂ may form a stem loop structure. In some embodiments, the primers may comprise two consecutive CCG triplets (CCG)₂ and additional nucleotides at the 5'- and 3'- end. In some embodiments, two sets of two consecutive CCG triplets (CCG)₂ are linked to each other through non-phosphodiester linkage. Exemplary non-nucleotide linkers include C4, C6, C8, C10 linkers, polyethylene glycol (PEG) linkers: PEG₄, PEG₆, PEG₁₀, PEG₁₂, peptide nucleic acid linkage, phosphorothoioester linkage. Exemplary arrangements of CCG triplets relative to each other are shown in Figure 15. Optionally the primers comprising one or more CCG triplets of Figure 15 may further comprise a sequence at the 5'-end that has at least 90% homology to the sequence directly upstream (i.e. at the 5'-end) of the CCG repeat tract. In some embodiments, such sequence can be at least 3, 5, 10 or more nucleotides.

[0108] In some embodiments, forward primers comprising one or more CCG triplets may further comprise an optional sequence at the 5'-end that has no homology with the target sequence of SEQ ID NO: 75. In some embodiments, the optional sequence at the 5'-end that has no homology to human genome. In some embodiments, forward primers comprising one or more CCG triplets may comprise a second portion and a third portion at the 5'-end of CCG triplets segment such that the second portion comprises a sequence that has at least 90% homology to the sequence directly upstream (i.e. at the 5'-end) of the CGG repeat tract and the third portion, which is at the 5'-end of the second portion, comprises a sequence that has no homology with the target sequence of SEQ ID NO: 75 (or any of the FMR2 gene).

[0109] Exemplary forward primer sequences comprising CCG triplets for detecting CCG repeat tract in the 5'-UTR of FMR2 gene can be any of the CCG triplet containing reverse primers used to detect CGG repeat tract in the 5'-UTR of FMR1 gene as listed in Table 2 above.

[0110] PCR Assay Condition:

[0111] Exemplary PCR assay mix may comprise 1X PCR buffer; 2mM MgCl₂; 6% DMSO (Sigma); 1.7% Q solution (Qiagen); 0.2 mM each of dGTP, dATP, dCTP, and dTTP; 1 unit of Roche FastStarTaq polymerase; and 0.6 μM each of the forward primer, reverse primer, and M13 reverse linker primer. The dNTP concentration can vary from 50-300 uM.

[0112] Optionally, dGTP may be partially or completely replaced with a dGTP analogue such as 7-deaza-2-deoxyGTP to improve the amplification efficiency of the highly GC rich

CGG repeat sequence. dGTP analogues are described in US patent 5,658,764 and 4,804,748. 7-deaza-dGTP prevents the Hoogsteen bond formation without interfering Watson-Crick base pairing and is incorporated by Taq polymerase (Seela et al. Biochemistry 1982; 21: 4338-4343). Several reagents may optionally be added to the PCR reaction mixture to improve amplification efficiency of the highly GC rich CGG repeat sequence of FMR1 gene and CCG repeat sequence of FMR2 gene. Exemplary reagents include but not limited to betaine, formamide, glycerol, dimethyl sulfoxide (DMSO). In one embodiment, the betaine concentration may be from 1M to 3M. In some embodiments, the betaine concentration can be 1.5M to 2.5M. In some embodiments, the DMSO concentration can be between 1%-3%. In some other embodiments, the DMSO concentration can be 1.5%. The PCR assay may further include controls such as positive controls and negative controls. Exemplary controls include nucleic acid samples with known status of CGG repeats for FMR1 gene and CCG repeats for FMR2 gene. Positive controls for normal or wild type FMR1 gene (i.e., less than 55 tandem CGG repeats), the premutation (55-200 tandem CGG repeats), and the full mutation (greater than 200 tandem CGG repeats) may be used. Similarly, for FMR2 gene, positive controls for normal or wild type FMR2 gene (i.e., less than 31 tandem CCG repeats), the premutation (60-200 tandem CCG repeats), and the full mutation (greater than 200 tandem CCG repeats) may be used.

[0113] Nucleic Acid Detection

[0114] Expansion of CGG repeats in the 5'-UTR of the FMR1 gene and CCG repeats in the 5'-UTR of the FMR2 gene can be detected by methods known in the art. Exemplary methods include but are not limited to sequencing, hybridization of specific probes such as Southern blot, detection by size such as by capillary electrophoresis, column chromatography or mass spectrometry. Detectable labels can be used to identify the probe hybridized to a nucleic acid such as genomic nucleic acid. Detectable labels include but are not limited to fluorophores, isotopes (*e.g.*, ³²P, ³³P, ³⁵S, ³H, ¹⁴C, ¹²⁵I, ¹³¹I), electron-dense reagents (*e.g.*, gold, silver), nanoparticles, enzymes commonly used in an ELISA (*e.g.*, horseradish peroxidase, betagalactosidase, luciferase, alkaline phosphatase), chemiluminiscent compound, colorimetric labels (*e.g.*, colloidal gold), magnetic labels (*e.g.*, DynabeadsTM), biotin, digoxigenin, haptens, proteins for which antisera or monoclonal antibodies are available, ligands, hormones, oligonucleotides capable of forming a complex with the corresponding oligonucleotide complement.

[0115] **Probes**

[0116] Oligonucleotide probes can be designed which are between about 10 and about 100 nucleotides in length and hybridize to the amplified region. Oligonucleotides probes are preferably 12 to 70 nucleotides; more preferably 15-60 nucleotides in length; and most preferably 15-25 nucleotides in length. The probe may be labeled. Amplified fragments may be detected using standard gel electrophoresis methods. For example, in preferred embodiments, amplified fractions are separated on an agarose gel and stained with ethidium bromide by methods known in the art to detect amplified fragments or may be detected by Southern blot.

[0117] Hybridization

[0118] The degree of complementarity between nucleic acid strands has significant effects on the efficiency and strength of hybridization between nucleic acid strands. This is of particular importance in detection methods that depend upon binding between nucleic acids. In one embodiment, the complementarity between the probe and the genomic nucleic acid may be "partial" in which only some of the nucleic acids' bases are matched according to the base pairing rules. In another embodiment, complementarity between the probe and the genomic nucleic acid may be "complete," "total," or "full".

[0119] The methods of the present invention can incorporate all known methods and means and variations thereof for carrying out DNA hybridization, see, e.g., Sambrook, et al., 1989, Molecular Cloning: A Laboratory Manual, Second Edition, Cold Spring Harbor Press, Plainview, NY.

[0120] The nucleic acid such as genomic nucleic acids and probes chosen for the nucleic acids are contacted under hybridization conditions. Hybridization conditions for nucleic acids in the methods of the present invention are well known in the art. For example, hybridization conditions may be high, moderate or low stringency conditions. Ideally, nucleic acids will hybridize only to complementary nucleic acids and will not hybridize to other non-complementary nucleic acids in the sample. The hybridization conditions can be varied to alter the degree of stringency in the hybridization and reduce background signals as is known in the art. For example, if the hybridization conditions are high stringency conditions, a nucleic acid will detectably bind to nucleic acid target sequences with a very high degree of complementarity. Low stringency hybridization conditions will allow for

hybridization of sequences with some degree of sequence divergence. The hybridization conditions will vary depending on the biological sample, and the type and sequence of nucleic acids. One skilled in the art will know how to optimize the hybridization conditions to practice the methods of the present invention.

[0121] Exemplary hybridization conditions are as follows. High stringency generally refers to conditions that permit hybridization of only those nucleic acid sequences that form stable hybrids in 0.018M NaCl at 65°C. High stringency conditions can be provided, for example, by hybridization in 50% formamide, 5X Denhardt's solution, 5X SSC (saline sodium citrate) 0.2% SDS (sodium dodecyl sulphate) at 42°C, followed by washing in 0.1X SSC, and 0.1% SDS at 65°C. Moderate stringency refers to conditions equivalent to hybridization in 50% formamide, 5X Denhardt's solution, 5X SSC, 0.2% SDS at 42°C, followed by washing in 0.2X SSC, 0.2% SDS, at 65°C. Low stringency refers to conditions equivalent to hybridization in 10% formamide, 5X Denhardt's solution, 6X SSC, 0.2% SDS, followed by washing in 1X SSC, 0.2% SDS, at 50°C.

[0122] Detectable label

[0123] The term "detectable label" as used herein refers to a molecule or a compound or a group of molecules or a group of compounds associated with a probe or primer and is used to identify the probe hybridized to a genomic nucleic acid or reference nucleic acid or the amplicon into which the primer incorporated.

[0124] Detectable labels include but are not limited to fluorophores, isotopes (e.g. 32P, 33P,35S, 3H, 14C, 125I, 131I), electron-dense reagents (e.g., gold, silver), nanoparticles, enzymes commonly used in an ELISA (e.g., horseradish peroxidase, beta-galactosidase, luciferase, alkaline phosphatase), chemiluminiscent compound, colorimetric labels (e.g., colloidal gold), magnetic labels (e.g., DynabeadsTM), biotin, digoxigenin, haptens, proteins for which antisera or monoclonal antibodies are available, ligands, hormones, oligonucleotides capable of forming a complex with the corresponding oligonucleotide complement.

[0125] In some embodiments, the detectable label is a fluorophore. Suitable fluorescent moieties include but are not limited to the following fluorophores working individually or in combination:

4-acetamido-4'-isothiocyanatostilbene-2,2'disulfonic acid; acridine and derivatives: acridine, acridine isothiocyanate; Alexa Fluor® 350, Alexa Fluor® 488, Alexa Fluor®

546, Alexa Fluor[®] 555, Alexa Fluor[®] 568, Alexa Fluor[®] 594, Alexa Fluor[®] 647 (Molecular Probes); 5-(2'-aminoethyl)aminonaphthalene-1-sulfonic acid (EDANS); 4-amino-N-[3vinylsulfonyl)phenyl]naphthalimide-3,5 disulfonate (Lucifer Yellow VS); N-(4-anilino-1naphthyl)maleimide; anthranilamide; Black Hole QuencherTM (BHQTM) dyes (biosearch Technologies); BODIPY dyes: BODIPY® R-6G, BOPIPY® 530/550, BODIPY® FL; Brilliant Yellow; coumarin and derivatives: coumarin, 7-amino-4-methylcoumarin (AMC, Coumarin 120),7-amino-4-trifluoromethylcouluarin (Coumarin 151); Cy2[®], Cy3[®], Cy3.5[®], Cy5[®], Cy5.5®; cyanosine; 4',6-diaminidino-2-phenylindole (DAPI); 5', 5"-dibromopyrogallolsulfonephthalein (Bromopyrogallol Red); 7-diethylamino-3-(4'-isothiocyanatophenyl)-4methylcoumarin; diethylenetriamine pentaacetate; 4,4'-diisothiocyanatodihydro-stilbene-2,2'disulfonic acid; 4,4'-diisothiocyanatostilbene-2,2'-disulfonic acid; 5-[dimethylamino]naphthalene-1-sulfonyl chloride (DNS, dansyl chloride); 4-(4'dimethylaminophenylazo)benzoic acid (DABCYL); 4-dimethylaminophenylazophenyl-4'isothiocyanate (DABITC); EclipseTM (Epoch Biosciences Inc.); eosin and derivatives: eosin, eosin isothiocyanate; erythrosin and derivatives: erythrosin B, erythrosin isothiocyanate; ethidium; fluorescein and derivatives: 5-carboxyfluorescein (FAM), 5-(4,6-dichlorotriazin-2yl) aminofluorescein (DTAF), 2',7'-dimethoxy-4'5'-dichloro-6-carboxyfluorescein (JOE), fluorescein, fluorescein isothiocyanate (FITC), hexachloro-6-carboxyfluorescein (HEX), QFITC (XRITC), tetrachlorofluorescein (TET); fluorescamine; IR144; IR1446; lanthamide phosphors; Malachite Green isothiocyanate; 4-methylumbelliferone; ortho cresolphthalein; nitrotyrosine; pararosaniline; Phenol Red; B-phycoerythrin, R-phycoerythrin; allophycocyanin; o-phthaldialdehyde; Oregon Green[®]; propidium iodide; pyrene and derivatives: pyrene, pyrene butyrate, succinimidyl 1-pyrene butyrate; QSY[®] 7; QSY[®] 9; OSY[®] 21; OSY[®] 35 (Molecular Probes); Reactive Red 4 (Cibacron[®] Brilliant Red 3B-A); rhodamine and derivatives: 6-carboxy-X-rhodamine (ROX), 6-carboxyrhodamine (R6G), lissamine rhodamine B sulfonyl chloride, rhodamine (Rhod), rhodamine B, rhodamine 123, rhodamine green, rhodamine X isothiocyanate, riboflavin, rosolic acid, sulforhodamine B, sulforhodamine 101, sulfonyl chloride derivative of sulforhodamine 101 (Texas Red); terbium chelate derivatives; N,N,N',N'-tetramethyl-6-carboxyrhodamine (TAMRA); tetramethyl rhodamine; tetramethyl rhodamine isothiocyanate (TRITC).

[0126] Other fluorescent nucleotide analogs can be used, see, e.g., Jameson et al. Meth. Enzymol. 1997; 278: 363-390; Zhu et al. Nucleic Acids Res. 1994; 22: 3418-3422. U.S. Patent Nos. 5,652,099 and 6,268,132 also describe nucleoside analogs for incorporation into

nucleic acids, e.g., DNA and/or RNA, or oligonucleotides, via either enzymatic or chemical synthesis to produce fluorescent oligonucleotides. U.S. Patent No. 5,135,717 describes phthalocyanine and tetrabenztriazaporphyrin reagents for use as fluorescent labels.

[0127] The detectable label can be incorporated into, associated with or conjugated to a nucleic acid. Label can be attached by spacer arms of various lengths to reduce potential steric hindrance or impact on other useful or desired properties. See, e.g., Mansfield et al. Mol. Cell. Probes. 1995; 9:145-156.

[0128] Detectable labels can be incorporated into nucleic acid probes by covalent or non-covalent means, e.g., by transcription, such as by random-primer labeling using Klenow polymerase, or nick translation, or, amplification, or equivalent as is known in the art. For example, a nucleotide base is conjugated to a detectable moiety, such as a fluorescent dye, e.g., Cy3TM or Cy5, TM and then incorporated into nucleic acid probes during nucleic acid synthesis or amplification. Nucleic acid probes can thereby be labeled when synthesized using Cy3TM- or Cy5TM-dCTP conjugates mixed with unlabeled dCTP.

[0129] Detection of nucleic acid by size:

[0130] Methods for detecting the presence or amount of polynucleotides are well known in the art and any of them can be used in the methods described herein so long as they are capable of separating individual polynucleotides by the difference in size of the amplicons. The separation technique used should permit resolution of nucleic acid as long as they differ from one another by at least one nucleotide. The separation can be performed under denaturing or under non-denaturing or native conditions--i.e., separation can be performed on single- or double-stranded nucleic acids. It is preferred that the separation and detection permits detection of length differences as small as one nucleotide. It is further preferred that the separation and detection can be done in a high-throughput format that permits real time or contemporaneous determination of amplicon abundance in a plurality of reaction aliquots taken during the cycling reaction. Useful methods for the separation and analysis of the amplified products include, but are not limited to, electrophoresis (e.g., agarose gel electrophoresis, capillary electrophoresis (CE)), chromatography (HPLC), and mass spectrometry.

[0131] In one embodiment, CE is a preferred separation means because it provides exceptional separation of the polynucleotides in the range of at least 10-3,000 base pairs with a resolution of a single base pair. CE can be performed by methods well known in the art, for

example, as disclosed in U.S. Pat. Nos. 6,217,731; 6,001,230; and 5,963,456, which are incorporated herein by reference. High-throughput CE apparatuses are available commercially, for example, the HTS9610 High throughput analysis system and SCE 9610 fully automated 96-capillary electrophoresis genetic analysis system from Spectrumedix Corporation (State College, Pa.); P/ACE 5000 series and CEQ series from Beckman Instruments Inc (Fullerton, Calif.); and ABI PRISM 3730 genetic analyzer (Applied Biosystems, Foster City, Calif.). Near the end of the CE column, in these devices the amplified DNA fragments pass a fluorescent detector which measures signals of fluorescent labels. These apparatuses provide automated high throughput for the detection of fluorescence-labeled PCR products.

- [0132] The employment of CE in the methods described herein permits higher productivity compared to conventional slab gel electrophoresis. By using a capillary gel, the separation speed is increased about 10 fold over conventional slab-gel systems.
- [0133] With CE, one can also analyze multiple samples at the same time, which is essential for high-throughput. This is achieved, for example, by employing multi-capillary systems. In some instances, the detection of fluorescence from DNA bases may be complicated by the scattering of light from the porous matrix and capillary walls. However, a confocal fluorescence scanner can be used to avoid problems due to light scattering (Quesada et al. Biotechniques. 1991; 10: 616-25.
- [0134] In some embodiments, nucleic acid may be analyzed and detected by size using agarose gel electrophoresis. Methods of performing agarose gel electrophoresis are well known in the art. See Sambrook et al., Molecular Cloning: A Laboratory Manual (2nd Ed.) (1989), Cold Spring Harbor Press, N.Y.
- [0135] In some embodiments, nucleic acid may be analyzed by HPLC. HPLC columns for separation of nucleic acid are available commercially such as HYDROCELL NS 1500, DEAE NP 10, QA NP 10 columns from Biochrom Labs Inc. (IN, USA) allow separation of nucleic acids. Ion pair reverse phase HPLC can also be used for DNA fragment analysis.

[0136] DNA sequencing:

[0137] In some embodiments, detection of nucleic acid is by DNA sequencing. Sequencing may be carried out by the dideoxy chain termination method of Sanger et al.(Proc. Natl. Acad.

Sci. USA 1977; 74: 5463-5467) with modifications by Zimmermann et al. Nucleic Acids Res. 1990; 18: 1067. Sequencing by dideoxy chain termination method can be performed using Thermo Sequenase (Amersham Pharmacia, Piscataway, NJ), Sequenase reagents from US Biochemicals or Sequatherm sequencing kit (Epicenter Technologies, Madison, Wis.). Sequencing may also be carried out by the "RR dRhodamine Terminator Cycle Sequencing Kit" from PE Applied Biosystems (product no. 403044, Weiterstadt, Germany), Taq DyeDeoxyTM Terminator Cycle Sequencing kit and method (Perkin-Elmer/Applied Biosystems) in two directions using an Applied Biosystems Model 373A DNA or in the presence of dye terminators CEQTM Dye Terminator Cycle Sequencing Kit, (Beckman 608000).

[0138] Alternatively, sequencing can be performed by a method known as Pyrosequencing (Pyrosequencing, Westborough, Mass.). Detailed protocols for Pyrosequencing can be found in: Alderborn et al. Genome Res. 2000; 10: 1249-1265.

[0139] Correlation of the number of CGG repeats in 5'-UTR of FMR1 gene and CCG repeats in 5'-UTR of FMR2 gene to Normal or Carrier/Affected Status

[0140] In some embodiments, individuals that are afflicted with a disease associated with an expansion of the CGG repeat region of the FMR1 gene can be distinguished from those that are normal. In some embodiments, a nucleic acid sample from the individual comprising the CGG repeat region is amplified using primers having one or more CCG triplets generate broad range of amplicons having varying lengths. The amplicons are separated into fractions according to size to have a distinctive pattern of tapering or "stutter". Male and female individuals can be differentiated into normal (5-44 CGG repeats), intermediate, premutation (55-200 CGG repeats), and full mutation alleles (> 200 CGG repeats) based on the size separation of the amplicons.

[0141] In some embodiments, individuals that are afflicted with a disease associated with an expansion of the CCG repeat region of the FMR2 gene can be distinguished from those that are normal. In some embodiments, a nucleic acid sample from the individual comprising the CCG repeat region is amplified using primers having one or more CGG triplets generate broad range of amplicons having varying lengths. The amplicons are separated into fractions according to size to have a distinctive pattern of tapering or "stutter". Male and female individuals can be differentiated into normal (6-30 CCG repeats), intermediate (31-60),

premutation (61-200 CCG repeats), and full mutation alleles (> 200 CCG repeats) based on the size separation of the amplicons.

[0142] Approaches to Patient Screening for FMR1 CGG Repeat Mutation Status and FMR2 CCG Repeat Mutation Status

In one embodiment, individuals may be tested for FMR1 CGG mutation status using the triplet primed amplification method of the invention in combination with other assays. In another embodiment, individuals may be tested for FMR2 CCG mutation status using the triplet primed amplification method of the invention in combination with other assays. One such assay also known as Capillary Southern Assay (CSA) was described by Strom et al. (Genet Med. 2007; 9(4):199-207) (see also US 20080124709 by Strom et al.). In the method described by Strom et al., gender of the nucleic acid sample is first confirmed by using an initial multiplex PCR that amplifies the FMR1 comprising the CGG triplet repeats, a region of the amelogenin gene that yields different fragment sizes from the X and Y chromosomes for gender determination, and a polymorphic region of the androgen insensitivity gene used as an internal control. The amplicons from the initial PCR reaction are then analyzed by an automated capillary DNA sequencer. After the initial PCR reaction and gender determination, the samples are further analyzed by the CSA method. The method involves restriction digestion of nucleic acid comprising the FMR1 gene and the 5'-UTR into multiple fragments and separating the fragments by capillary electrophoresis. The presence or absence of FMR1 gene is detected in the size separated fragments by gene specific PCR. The size of the fractions is indicative of CGG repeat expansion.

[0144] The triplet primed amplification method of the invention is useful for identifying larger premutations and full mutations. The CSA method enables accurate sizing of repeat number in the normal, intermediate, and small premutation zone as well as gender identification/confirmation. The combination of both the CSA method and triplet primed amplification for each individual would yield results for both males and females, and this would serve the purpose for carrier and/or newborn screening. The triplet primed amplification method should significantly reduce the number of Southern blot analyses to be performed as the method clearly distinguishes normal homozygous females, which represent ~95% of all female samples from full mutation carriers. A schematics of the two different assays run simultaneously for determining CGG repeat expansion is shown in Figure 17.

[0145] The CSA method of Strom et al. can also be applied to the detection and confirmation of CCG repeat mutation status of FMR2 gene. In some embodiments, the FMR1 CGG mutation status and the FMR2 CCG mutation status can be determined simultaneously in a single assay. In another embodiment, the FMR1 CGG mutation status and the FMR2 CCG mutation status can be determined using the same sample using two separate assays.

[0146] Estimation of premutation and full mutation frequency in FMR1 gene

[0147] In some embodiments, estimation of the number of CGG triplets in the CGG repeat tract of FMR 1 gene will allow an artisan to determine the premuation carrier frequency and full mutations in a population. The estimation of the number of CGG triplets in the CGG repeat tract of FMR 1 gene will also allow to calculate the incidence of FXS, FXTAS and FXPOI in a population. Premuation carrier frequency in females can be estimated from the experimental results and can be expressed as the fraction of female population tested to have premutation (≥55 CGG repeats) in FMR1 gene. The premuation carrier frequency in males as well as the full mutation frequency in both males and females can be calculated from the premuation carrier frequency in females. Full mutation frequency in an individual (male or female) can be calculated based on the premutation carrier frequency of a female using the following equation:

Full mutation frequency = premutation carrier frequency in females x S x 0.5

where S is the rate of expansion of a premutation into a full mutation and was calculated to be 0.107 (i.e. $\sim 10.7\%$ of premutations expand to full mutations), while 0.5 is the probability of inheriting either allele.

[0148] Premutation carrier frequency in males can be calculated using the equation:

Premut. carrier frequency in males = premut. carrier frequency in females x (1-S) x 0.5

[0149] Estimation of the incidence of FXTAS and FXPOI in a populations

[0150] In some embodiments, the incidence of FXTAS and FXPOI in a populations can be calculated by utilizing the estimated frequency of premutations in males or females, the reported penetrance of alleles in male and female carriers, and frequency of larger premutation alleles. The incidence of FXTAS can be calculated using the following equation:

incidence of FXTAS = premutation carrier frequency in males x P_{FXTAS} x frequency of alleles ≥ 70 CCG repeats

where P_{EXTAS} is penetrance in patients of larger premutation alleles (\geq 70 CGG repeats) in FXTAS estimated to be 1/3 (0.33).

[0151] The incidence of FXPOI incidence can be calculated using the following equation:

incidence of FXPOI = premutation carrier frequency in females x P_{FXPOI} x frequency of alleles ≥ 70 CCG repeats

where P_{FXPOI} is the penetrance of large premutation alleles (\geq 70 CGG repeats) in FXPOI estimated to be \sim 20%-28%.

EXAMPLE 1

Samples

[0152] Residual archived DNA from patient samples already processed for fragile X analysis using existing assays were de-identified and used in the development of the PCR method. DNA was extracted using the Qiagen M96 robot and reagents (Qiagen, Carlsbad, CA), Qiagen Gentra robot and reagents, and Qiagen 9604 reagents. The samples included a homozygous normal female (28/28 repeats), a heterozygous normal female (18/21), a premutation female (31/64), full mutation females (19/450 and 29/530), a normal male (28), a premutation male (74), and mosaic full mutation male (105/350). Once conditions for the PCR were established, a larger set of 47 anonymized, previously genotyped samples (Table 4), weighted toward full mutation carriers (male or female), were used to test the robustness of the assay.

Table 4. Panel of Previously *FMR1* Genotyped Samples Used to Test Performance of the Triplet-Primed PCR/Capillary Electrophoresis-Based Assay.

Sample #	Previous Genotype Call	Southern Blot Analysis	Classification by Triplet-Primed PCR/CE
1	23/800		FF
2	29/590		FF
3	687		FM
4	19/88		PF
5	104		PM

Sample #	Previous Genotype Call Southern Blot Analysis		Classification by Triplet-Primed PCR/CE
6	82		PM
7	28/400		FF
8	506-1006		FM
9	29/286		FF
10	28/98		PF
11	660		FM
12	400		FM
13	29/29		NF
14	696		FM
15	67/500		FM
16	28/29		NF
17	78/2000	78 (60%, unmethylated) / 2000 (40%, methylated)	FM
18	45		GM
19	QC Blank		No Amplification
20	23/400-1250		FF
21	650		FM
22	850		FM
23	45/600		FF
24	30/550		FF
25	30/49		GF
26	713		FM
27	500		FM
28	400		FM
29	30/520		FF
30	19/19		HF
31	500		FM
32	30/30		HF
33	32/73		PF
34	31/480		FF
35	170-240		FM
36	400		FM
37	600		FM
38	200-500	200 (10%, unmethlyated) / 200-500 (90%, methylated)	FM

Sample #	Previous Genotype Call	Southern Blot Analysis	Classification by Triplet-Primed PCR/CE
39	420		FM
40	30/74		PF
41	853		FM
42	30/313		FF
43	606		FM
44	433		FM
45	140/(200- 400)	140 (5-10%, unmethylated) / 200-400 (90-95%, methylated)	FF
46	300-500		FM
47	590-1143		FM
48	666		FM

NF: normal female; GF: intermediate female; PF: premutation female; FF: affected female; NM: normal male; GM: intermediate male; PM: premutation male; FM: affected male.

[0153] To evaluate assay sensitivity, three DNA sets were made using archived, previously genotyped samples, each set harboring two DNA samples mixed at different ratios; set # A contained DNA samples from two males, one from a normal male with 28 CGG repeats, and one from affected male with 850 CGG repeats; set B contained two males, a normal male (28 CGG) and affected male (853), and set C contained two female samples, a normal (28/28) and an affected female (19/450).

[0154] To determine whether assay methods could detect fragile X CGG mutation status on blood spots, residual blood from previously genotyped patients were spotted on Guthrie cards. Blood spots were punched from each card using BSD1000 GenePunch Instrument to generate two punches (3.2 mm each) per patient and DNA was extracted using Qiagen BioSprint reagents. DNA from blood spots was analyzed by the triplet-primed PCR/capillary electrophoresis as described in Example 2. A total of 37 samples containing the different *FMR1* CGG expansions in males and females were tested (Table 5). The series contained fourteen samples with intermediate repeats, nine samples with pre-mutations, and four samples with full mutations (as determined by Southern blot analysis), one female (20/500 repeats) and three males (300-700; 1250; 1250 repeats).

Table 5. Samples used for blood spots.

Sample #	Genotype*	Previous	Classification by Triplet-
		Result	Primed PCR/CE
1	29/46	IF	No Expansion
2	23/55	PF	PF
3	28/66	PF	PF
4	31/53	IF	No Expansion
5	20/500	FF	FF
6	60	PM	PM
7	46	IM	No Expansion
8	53	IM	No Expansion
9	28/28	HF	No Expansion
10**	-	-	No Amplification
11	31/45	IF	No Expansion
12	30/49	IF	No Expansion
13	29/29	HF	No Expansion
14	19/54	IF	No Expansion
15	30	NM	No Expansion
16	34/50	IF	No Expansion
17	300-700	FM	FM
18	44/57	PF	PF
19	32/45	IF	No Expansion
20**	-	-	No Amplification
21	50	IM	No Expansion
22	20/51	IE	No Europaise
22	30/51	IF	No Expansion
23	29/70	PF	PF
24	1250	FM	FM

25	46	IM	No Expansion
26	30/30	HF	No Expansion
27	29/78	PF	PF
28	29/29	HF	No Expansion
29	30/64	PF	PF
30**	-	-	No Amplification
31	1250	FM	FM
32	35/59	PF	PF
33	30/30	HF	No Expansion
34	30/66	PF	PF
35	30/30	HF	No Expansion
36	29/29	HF	No Expansion
37	30/30	HF	No Expansion
38	30/30	HF	No Expansion
39	32/47	IF	No Expansion
40**	-	-	No Amplification
41	30/49	IF	No Expansion

Genotype is listed as the number of triplet repeats on each *FMR1* allele. Abbreviations are as in Table 4. No Expansion includes normal and intermediate alleles. ** Samples 10, 20, 30 and 40 contained blank cards for cross contamination detection. These samples showed no amplification (NA).

[0155] To determine the premutation frequency in US population, residual anonymized DNA samples submitted for Cystic Fibrosis (CF) carrier screening and samples submitted for screening for diseases more prevalent in Ashkenazi Jews were used. A total of 13770

samples were analyzed in this study, >99% of the sample were from females. Table 6 shows the breakdown of samples by self-identified ethnicity. Greater than 43% did not provide ethnicity information, while 4.9% of samples listed as other or mixed ethnicities. Caucasian accounted for 20.8% of the total samples, Hispanics accounted for 12.5% of samples, samples from each of African Americans and Asians accounted for 1.8% of the total samples.

Table	5: '	Composition	of	the	popu.	lation	studied.	
-------	------	-------------	----	-----	-------	--------	----------	--

Ethnicity	Samples Screened	% Patients
Caucasian	2866	20.8
African American/Black	247	1.8
Hispanic	1723	12.5
Ashkenazi Jewish	2011	14.6
Asian	254	1.8
Other	674	4.9
Not given	5995	43.5
Total	13770	100.0

EXAMPLE 2 Assay Development and Assay Conditions

[0156] The CGG repeats in the FMR1 gene were amplified by PCR. The PCR reaction included a fluorescently labeled forward primer and reverse primers comprising at least two CCG triplets such that the reverse primers can anneal to CGG repeat sequence. In one example, the reverse primer also comprised the junction sequence of the 3'-end of the CGG repeat tract and flanking sequence and a portion of the sequence downstream of CGG repeat. The reverse primers having at least two CCG triplets and without such junction sequence will hybridize randomly inside the CGG repeat sequence and initiate DNA synthesis. This will result in a range of amplicons of different sizes (stutter). The reverse primers further comprised a M13 linker sequence at its 5'-end such that a M13 primer can anneal to this linker sequence and further amplify the amplicons thus generated.

[0157] The PCR mix for amplification of all normal, premutation, and full mutations of the CGG repeats were optimized by formulating various combinations of mixes using Roche Expand polymerase PCR mix buffer 2 or 3, Roche High-GC PCR mix, and incorporating

various additives at different concentrations including Roche Q-solution, DMSO, HiDi Formamide, Tween-20, or combinations thereof. Initially, Roche Expand kit without the addition of deaza-dGTP was tested. However, Roche FastStarTaq kit reagents (Roche Applied Science, Indianapolis, IN) provided better performance in terms of signal intensity under identical cycling conditions.

[0158] The PCR mix contained 1X PCR buffer; 2mM MgCl₂; 6% DMSO (Sigma); 1.7% Q solution (Qiagen); 0.2 mM each of 7-deaza-2-deoxyGTP, dATP, dCTP, and dTTP; 1 unit of Roche FastStarTaq polymerase; and 0.6 μM each of the forward primer (FMR1F), reverse primer, and M13 reverse linker primer. Several variations of the reverse primers were tested together with the FMR1F forward primer. The sequences of the primers are listed below in Table 7.

Table 7: Primers used in the PCR reaction

Primer	Sequence	SEQ ID NO:
FMR1F Forward	6-FAM_TGTAAAACGACGGCCAGTGCTCAGCTCCGTTTCGGTTTCA CTTCCGGT	2
FMR1R Reverse	<u>CAGGAAACAGCTATGACC</u> ctcgaggcccaGCCGCCGCCGCC	3
FMR1CCGR Reverse	<u>CAGGAAACAGCTATGACC</u> CCGCCGCCGCC	4
CGG1 Reverse	<u>CAGGAAACAGCTATGACC</u> CGCCGCCCC	5
CGG2 Reverse	<u>CAGGAAACAGCTATGACC</u> CCGCCGDCGCCGCCC	6
CGG3 Reverse	<u>CAGGAAACAGCTATGACC</u> CCGCCGCCGCCGCCGCCGCCGCCGCCG	7
CGG4 Reverse	<u>CAGGAAACAGCTATGACC</u> CCGCCGCCGCCGCCGCCGCCG	8
CGG3-1 Reverse	<u>CAGGAAACAGCTATGACC</u> CCGCCGCCGCCGCCGCCGCCGCCGCC	9
CGG4-1 Reverse	<u>CAGGAAACAGCTATGACC</u> CCGCCGCCGCCGCCGCC	10
FMR1R_7 Reverse	<u>CAGGAAACAGCTATGACC</u> CTCGAGGCCCAGCCGCCGCCGCC	71
FMR1R_8 Reverse	<u>CAGGAAACAGCTATGACC</u> CTCGAGGCCCAGCCGCCGCCGCCGCC	72
FMR1R_9 Reverse	<u>CAGGAAACAGCTATGACC</u> CTCGAGGCCCAGCCGCCC	73
FMR1R_10 Reverse	<u>CAGGAAACAGCTATGACC</u> CTCGAGGCCCAGCCGCCG	74
M13 reverse linker primer	CAGGAAACAGCTATGACC	11

D: A, G or T

[0159] Three microliters of genomic DNA (10-150 ng total) was added for a final PCR reaction volume of 15 μ l. The cycling program, performed on an ABI 9700 thermal cycler,

was started by incubating the mix at 98°C for 10 minutes to activate the polymerase, followed by 10 cycles of denaturation at 97°C for 35 sec, annealing at 64°C for 2 minutes, and extension at 68°C for 8 minutes, followed by 25 cycles of denaturation at 97°C for 35 sec, annealing at 64°C for 2 minutes, and extension at 68 °C for 8 min 20 sec (with 20 sec extension each additional cycle). For all 35 cycles, the ramp rate was adjusted to 64% for the denaturation and extension steps (ramping up at 0.5°C), and to 25% for the annealing step (ramping down at 0.5°C). Two microliters of the PCR product (some times diluted 1:5 in H₂O) were mixed with HiDi Formamide and Map marker-1000 size standard (Bioventures Inc., Murfreesboro, TN) and samples were injected on an ABI 3730 DNA Analyzer equipped with a 36-cm capillary loaded with Pop-7 polymer (Applied Biosystems) at 3 volts for 7 seconds. Data were analyzed using GeneMapper software (Applied Biosystems).

[0160] The forward primer FMR1F is fluorescently labeled with FAM and anneals to 5'-UTR of FMR1 gene upstream of the CGG repeats. The reverse primers were designed to comprise at least two consecutive CCG triplets such that the primers anneal to the CGG repeat sequences of the FMR1 gene. Reverse primer FMR1R (SEQ ID NO: 3) comprise 4 CCG triplet sequences and the junction sequence directly next to the 3'-end of the CGG repeat tract, the later indicated in bold lowercase. The reverse primers also comprise M13 linker sequence at its 5'-end (underlined) which provides hybridization sites for M13 reverse linker primer in further round of amplification. Reverse primers FMR1CCGR, CGG1, CGG2, CGG3, CGG4, CGG3-1, CGG4-1, FMR1_7, FMR1_8, FMR1_9, FMR1_10 can hybridize randomly to the CGG repeat sequence of FMR1 gene and initiate DNA synthesis to generate a range of amplicons of different sizes (stutter). A schematics of the PCR assay indicating the forward and reverse primers binding sites relative to the CGG repeat sequences is shown in Figure 2.

[0161] The inclusion of the M13 linker in the reverse primers also helped to boost the signal by about 30%. The incorporation of ramping in the cycling program helped further boost signal amplification by at least three fold.

EXAMPLE 3

Effect of the addition of 7-deaza-2-deoxyGTP on PCR Amplification of CGG Repeats

[0162] The effect of substituting dGTP with 7-deaza-2-deoxyGTP in the PCR reaction was tested to improve the amplification efficiency of the highly GC rich CGG repeat sequence. Deaza-dGTP prevents the Hoogsteen bond formation without interfering Watson-Crick base pairing and is incorporated by Taq polymerase (Seela et al. Biochemistry 1982; 21: 4338-4343). The effect of increasing concentration of 7-deaza-2-deoxyGTP was tested in the PCR reaction mixture using genomic DNA isolated from a heterozygote female carrier (29/530). Increasing the concentration of 7-deaza-2-deoxyGTP not only boosted the signal from stutter amplification, but also led to the appearance of a peak that migrated very late on the ABI 3730 capillary and was evident only when a full mutation was present (Fig 3 and Fig 4). Even though amplicons containing 7-deaza-2-deoxyGTP do not stain with ethidium bromide, the presence of a premutation and a full mutation could be detected in a single PCR reaction by stutter and full mutation signals, thus distinguishing homozygous normal females from heterozygous full mutation-carrier females.

[0163] Using a panel of eight anonymized, previously genotyped samples, various possibilities of CGG triplet repeat status were tested. The assay could distinguish between normal and intermediate alleles (< 55 CGG repeats) from premutation and full mutation alleles (with ≥ 55 CGG repeats) as shown in the electropherogram (Figure 4). Normal females and normal males display alleles only in the gray (or left hand shaded) area of the electropherogram, while premutation males and females, display alleles in the pink area (or right hand shaded), with stutter amplification prior to the actual premutation peak. Full mutation carrying males and females display, in addition to normal allele in females, a tapering stutter amplification followed by a signal at the end of electropherogram that is present only whenever a full mutation is present. Most significantly, normal homozygote females can be distinguished from full mutation carrier females using this method (Figure 4).

EXAMPLE 4

Amplification of CGG Repeats by Different Reverse Primers

[0164] The efficiency of amplification of CGG repeats by the reverse primers were tested using a PCR reaction mixture comprising 1X PCR buffer; 2mM MgCl₂; 6% DMSO (Sigma); 1.7X Q solution (Qiagen); 0.2 mM each of 7-deaza-2-deoxyGTP, dATP, dCTP, and dTTP; 1

unit of Roche FastStarTaq polymerase; and 0.6 μM each of the forward primer (FMR1F), M13 reverse linker primer, and one of the reverse primers: FMR1R, CGG1, CGG2, CGG3, CGG4, CGG3-1 and CGG4-1 using genomic DNA isolated from 4 samples: full mutation female (FF 23/800) full mutation male (FM 687), premutation female (PF 19/88) and homozygous female (HF 29/29).

[0165] All reverse primers generated a stutter pattern of amplicons for the three samples: full mutation female, full mutation male and premutation female (Figures 8-14). Reverse primers FMR1R, CGG1, CGG2, CGG4-1 generated full mutation peaks which appeared around the 1050 bases for samples FF 23/800 and FM 687 as indicated in the electropherogram of figures 8, 9, 10 and 14.

EXAMPLE 5

Effect of the Addition of Various Concentrations of GTP, 2-deoxyGTP and 7-deaza-2-deoxyGTP on PCR Amplification of CGG Repeats by Different Reverse Primers at Different Annealing Temperatures

The effect of increasing concentration of GTP, 2-deoxyGTP, and 7-deaza-2-deoxy [0166]GTP on the amplification efficiency of the highly GC rich CGG repeat sequence were tested. Varying amounts of GTP, 2-deoxyGTP, and 7-deaza-2-deoxyGTP were added to the PCR reaction mixture comprising 1X PCR buffer; 2mM MgCl₂; 6% DMSO (Sigma); 1.7X Q solution (Qiagen); 0.2 mM each of dATP, dCTP, and dTTP; 1 unit of Roche FastStarTaq polymerase; and 0.6 µM each of the forward primer (FMR1F), M13 reverse linker primer, and one of the reverse primers: FMR1R, FMR1R 7, FMR1R 8, FMR1R 9, and FMR1R 10 using genomic DNA isolated from a heterozygote female carrier at different annealing temperatures. The amplicons were separated into fractions according to size to have a distinctive pattern of tapering or "stutter". Full mutation alleles (> 200 CGG repeats) were determined based on the size separation of the amplicons. A reference assay condition was arbitrarily determined to comprise 1X PCR buffer; 2mM MgCl₂; 6% DMSO (Sigma); 1.7X Q solution (Qiagen); 0.2 mM each of 7-deaza-2-deoxyGTP, dATP, dCTP, and dTTP; 1 unit of Roche FastStarTaq polymerase; and 0.6 µM each of the forward primer (FMR1F), M13 reverse linker primer, and reverse primer FMR1R at an annealing temperature of 64°C. Results from each of the amplification conditions were compared to the reference condition to determine the efficiency of detecting stutter pattern and full mutation.

[0167] The addition of GTP to the reaction mixture up to a final concentration of 200μM did not have a significant effect on the amplification efficiency of the highly GC rich CGG repeat sequence. The amplification efficiencies were comparable or better in reaction conditions having 7-deaza-2-deoxyGTP concentration of .85 mM to 1 mM. The amplification efficiencies with reverse primer FMR1R_9 under various reaction conditions were comparable or better than the results obtained from the reference condition, while the amplification efficiencies with reverse primer FMR1R_8 under various reaction conditions tested were worse than the results obtained from the reference condition. A comparison of the results of the different reaction conditions to the reference condition is presented in Table 8 below.

Table 8: Comparison of amplification efficiency for detection of stutter pattern and full mutation in presence of GTP, 7-deaza-2-deoxyGTP, using various reverse primers

	Primer	Primer	Length	Deaza-	dGTP	СТР	0/0	Annealing			
Condition	Forward		of	dGTP			⁷⁰ dGTP/GTP	Temp	Result	Stutter	Full Mutation
	roiwaiu	Reverse	CCG	uM	uivi	uivi	uG1P/G1P	(°C)			
Ref	FMR1F	FMR1R	3 2/3	200	0		0.0	64	Reference	Present	Present
Ref	I WIKII	I WIKIK	2/3	200			0.0		Condition	resent	i resent
1	FMR1F	FMR1R 9	2 2/3	650	12		1.8	60	Comparable	Present	Present/Strong
	1111111		2 2, 3		12		1.0		or better	resent	recent Suong
2	FMR1F	FMR1R_9	2 2/3	650	12		1.8	64	Better	Present	Present/Strong
3	FMR1F	FMR1R 9	2 2/3	1000	12		1.2	L 60	Comparable	Present	Present/Strong
		, ,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,		1000	12		1.2		or better	resent	resem suong
4	FMR1F	FMR1R_9	2 2/3	1000	12		1.2	64	Better	Present	Present/Strong
5	FMR1F	FMR1R_9	2 2/3	850	6		0.7	60	Comparable	Present	Present/Strong
				000			5.,		or better	100011	2 10001111 2 11 0 11 9
6	FMR1F	FMR1R_9	2 2/3	850	6		0.7	64	Better than	Present	Present/Strong
)							Std		
7	FMR1F	FMR1R 9	2 2/3	850	2		0.2	60	Comparable	Present	Present/Strong
		_							or better		
		FMR1R_9	2 2/3	850	2		0.2	64			Present/Strong
9	FMR1F	FMR1R	3 2/3	500	12		2.3	64	Worse	Present	Absent/week
10	FMR1F	FMR1R	3 2/3	650	12		1.8	64	Worse	Present	weak
11	FMR1F	FMR1R	3 2/3	1000	12		1.2	64	Comparable	Present	Present
12	FMR1F	FMR1R	3 2/3	200		12	5.7	64	Comparable	Present	Present
13	FMR1F	FMR1R_10	2	500		12	2.3	60			
14	FMR1F	FMR1R	3 2/3	500		12	2.3	64			

15	FMR1F	FMR1R	3 2/3	850	6	0.7	64	Comparable or better	Present	Present
16	FMR1F	FMR1R	3 2/3	850	2	0.2	64	Comparable or better	Present	Present
17	FMR1F	FMR1R_7	4 2/3	650	12	1.8	70	Comparable or better	Present	Present
17	FMR1F	FMR1R_7	4 2/3	1000	12	1.2	67	Comparable or better	Present	Present
19	FMR1F	FMR1R_7	4 2/3	200	0	0.0	64	Worse	Present	Absent/weak
20	FMR1F	FMR1R_8	5 2/3	500	12	2.3	64	Worse	Present	Absent/weak
21	FMR1F	FMR1R_8	5 2/3	200	0	0.0	64	Worse	Present	Absent/weak

EXAMPLE 6Assessment of the Robustness of the Assay to Detect CGG Repeats

[0168] To assess the robustness of the PCR assay of Example 2, previously genotyped DNA samples were tested using primers FMR1F and FMR1R, M13R in presence of 7-deazadGTP. The person performing the test was blinded to the status of CGG repeats in this panel, and the results scored by a second person. The panel of anonymized forty-seven samples contained homozygous females (N=3), normal heterozygous females (N=1), intermediate female (N=1), premutation females (N=4), and full mutation females (N=11), in addition to intermediate male (N=1), premutation male (N=2), and full mutation males (N=24) (Table 2) The panel assembled was heavy on full mutation carrier males and females. All the samples were identified correctly showing the capability of the new method in identifying presence of full mutations, and distinguishing them from normal homozygous females, or normal heterozygous females with two close alleles (samples H2, 28/29).

EXAMPLE 7

Detection of Mosaics

[0169] The ability to detect in a patient blood sample mosaic alleles of premutation and full mutations identified were tested using the PCR assay discussed in Example 2. Figure 4 shows the results of a mosaic male with a premutation allele of 105 repeats and a full mutation allele of 350, identified from a blood sample. Mosaic samples listed in Table 2 with mosaic full mutations including two male samples (samples 15, 17) and two female samples (samples 2, 7) were tested. Sample 17, which was mosaic for 78 and 2000 repeats showed amplification signal for the longest full mutation (Figure 5). The signal for the full mutation

was weak but detectable. Previous Southern blot analysis of this sample showed the unmethylated premutation allele accounted for 60% of alleles while the methylated full mutation 2000 repeat allele accounted for the rest in this affected male sample. The method also detected mosaic premutations in full mutation carrier females (Figure 5). These results show that the new PCR method is able to detect mosaic alleles in both males and females.

EXAMPLE 8

Assay Sensitivity

[0170] Assay sensitivity was tested by mixing DNA samples from males and females with normal of full mutations creating artificial mosaics. Three DNA sets were prepared using archived, previously genotyped samples, each set harboring two DNA samples mixed at different ratios; set "A" contained DNA samples from two males, one from a normal male with 28 CGG repeats, and one from affected male with 850 CGG repeats; set "B" contained two different male DNAs, a normal male (28 CGG repeats) and affected male (853 CGG repeats), and set "C" contained two female samples, a normal (28/28 repeats) and an affected female (19/450 repeats). Figure 6 shows the results from detecting these samples. For male DNA, the full mutation was detectable at levels of 5% in set "A" and 10% in set "C". As shown in Figure 6, the intensity of the stutter amplification and the full mutation signals increase as the level of full mutation increases. In set "B", from mixed male samples, it is interesting to note the ability of the assay to detect as low as 11% mosaic full mutation content in female DNA that contained alleles of three lengths (19, 28, and 450 repeats), where the smaller, normal alleles would have been favored during PCR. At 10% mosaic content, a weak signal from a full mutation is present (Fig 6C), while stutter amplification is appearing reliably at 20% content of the full mutation. On average, full mutation expansion could be detected by both the stutter amplification and full mutation signal detection, at levels as low as 5% to 10% in both males and females.

EXAMPLE 9

Concordance Study

[0171] The assay results of new PCR method were compared against the existing assay method which is composed of a two step-process: a PCR incorporating gender marker, and a subsequent Southern blot analysis for any sample with an expansion length of 45 repeats or above, as well as any apparent homozygous females. Previously genotyped 1,275 whole blood samples were analyzed using the new PCR mix including PCR primers FMR1F and FMR1R, M13R and 7-deaza dGTP; the results indicated 100% concordance with the existing

method. As shown in Table 9, of these 1,275 samples a total of six patient samples (three female, three male) harbored full mutations.

Table 9. Concordance of Triplet-Primed PCR/Capillary Electrophoresis Method with Existing PCR/Southern Blot Method for Determination of *FMR1* Trinucleotide Expansion Status

Classification	Number
Normal (HF, HetF, NM)	1237
Intermediate Female	15
Intermediate Male	7
Premutation Female	10
Premutation Male	0
Full Mutation Female	3
Full Mutation Male	3

HF, normal homozygous female; HetF, normal heterozygous female; NM, normal male.

EXAMPLE 10

Detection From Blood Spots and Saliva

[0172] The CGG repeat-primed PCR assay was tested for the detection of CGG repeats in FMR1 gene from blood spots, since newborn screening programs rely on blood spotted on Guthrie cards. Residual blood from previously genotyped patients were spotted on Guthrie cards. Blood spots were punched from each card using BSD1000 GenePunch Instrument to generate two punches (3.2 mm each) per patient and DNA was extracted using Qiagen BioSprint reagents. DNA from blood spots was analyzed by the new PCR method. A total of 37 samples containing the different *FMR1* CGG expansions in males and females were tested (Table 5). Of the 37 samples, fourteen samples had intermediate repeats, nine samples had pre-mutations, and four samples had full mutations (as determined by Southern blot analysis): one female (20/500 repeats) and three males (300-700; 1250; 1250). Figure 7 shows an example of the results obtained. All the samples where genotyped correctly by the new qualitative method as harboring normal, intermediate, pre-mutation or full mutations.

EXAMPLE 11

Estimation of Premutation Carrier Frequency in US Population

A total of 13770 genomic DNA samples were analyzed by CCG repeat-primed PCR assay for the presence of premutations and/or full mutation. The DNA samples were obtained from individual samples submitted for CF screening and screening diseases prevalent among Ashkenazi Jews. The PCR reaction mixture comprised 1X reaction buffer, MgCl₂ at final concentration of 2 mM, 1 unit of FastStarTaq polymerase (all from Roche Applied Science, Indianapolis, IN), 0.6 µM each of FMR1F, FMR1R and M13R primers, 0.2 μM each of the nucleotides dATP, dCTP, dTTP and 7-deaza-2-deoxy GTP, 1.7X Q solution (Qiagen, Carlsbad, CA) and 6% DMSO (Sigma, St. Louis, MO). The PCR reaction products were separated on ABI 3730 DNA Analyzer (Applied Biosystems, Carlsbad, CA) and analyzed using GeneMapper Software (Applied Biosystems, Carlsbad, CA) as described above. Samples showing fragments containing ~50 or more CGG repeats were further analyzed using the Capillary Southern Assay (CSA) was described by Strom et al. (Genet Med. 2007; 9(4):199-207) (see also US 20080124709 by Strom et al.). The method includes a multiplex PCR assay together with capillary electrophoresis (CE) and Southern blot to determine the gender of the individual, the size of the CE fractions which is indicative of CGG repeat expansion.

- [0174] Figure 18A provides an example of the results obtained for a CCG repeat-primed PCR assay. The sample was identified as harboring a non expanded allele (< 55 repeats, indicated) and an expanded allele (> CGG 55 repeats, indicated). The sample was reanalyzed using CSA method with gender confirmation. As shown in Fig 18B, the sample was confirmed to be from a female as only the X-chromosome copy of the Amelogenin gene was amplified. The sample showed heterozygosity for 29 and 65 repeats.
- [0175] A total of 11759 samples from CF-screened females were analyzed as a representative of the US population. Forty eight premutations (> 55 CGG repeats) were identified from the CF-screened samples showing a premutation carrier frequency of 1:245. The distribution of the 48 premutation alleles identified is shown in Fig 19. Twelve of the 48 alleles (25%) harbored ≥70 CGG repeats.
- [0176] One female sample from a Caucasian patient showed evidence of full mutation as shown by stutter pattern of amplification and detection of the late migrating fragment on capillary electrophoresis (Fig 20A). That sample was further analyzed by Capillary Southern

Assay and was shown to harbor 30/(400-800) repeats, confirming the finding of CCG triplet primed PCR. The methylation status of the DNA sample was analyzed by digestion with EcoRI and XhoI restriction endonuclease followed by agarose gel electrophoresis. The sample showed the appearance of additional DNA fragments (indicated by arrow) along with the expected sizes indicating heterogeneity in DNA methylation (Fig 20B).

[0177] The distribution of premutation carrier frequency among different US sub population or ethnic groups were estimated by analyzing a total of 13770 samples including DNA samples submitted for Cystic Fibrosis (CF) carrier screening and samples submitted for screening for diseases more prevalent in Ashkenazi Jews. Premutation carrier frequency in females was estimated from the experimental results of the CCG primed PCR assay and the premutation carrier frequencies in males were calculated from the estimated premuation carrier frequency in females.

Premutation carrier frequency in males can be calculated using the equation:

Premutation carrier frequency in males = premutation carrier frequency in females X (1-S) X 0.5

where S is the rate of expansion of a premutation into a full mutation and was calculated to be 0.107 (i.e. $\sim 10.7\%$ of premutations expand to full mutations), while 0.5 is the probability of inheriting either allele.

[0178] A total of 63 premutations were identified among the 13770 samples resulting in a premutation carrier frequency of 1:219. This frequency is higher due the over representation of the samples from Ashkenazi Jews. Therefore, taking the CF-screened population as approximate representative of the US population, the premutation carrier frequency was found to be 1:245 (48 premutation samples identified out of 11759 CF-screened females). Table 10 shows the results from the various self-identified ethnic groups within CF population samples tested. Premutation carrier frequency among 2866 Caucasians was 1:169, while it was 1:287 in 1723 Hispanics. No premutations were identified in Asians, while the frequency in African American was 1:124. In Ashkenazi Jewish samples, 15 premutations were detected from 2011 DNA samples with a premutation carrier frequency of 1:134. Only one tested sample from Ashkenazi Jews harbored a premutation ≥70 CGG repeats (80 repeats).

Table 10: Distribution of premutation carrier frequency in different US subpopulation

Ethnicity	# Samples Screened	# Premutations Identified	Frequency	
Caucasian	2866	17	1:169	
Hispanic	1723	6	1:287	
African American/Black	247	2	1:124	
Ashkenazi Jewish	2011	15	1:134	
Asian	254	0	N/A	
Other	674	2	1:337	
Not Given	5995	21	1:285	
Total	13770	63	1:219	

[0179] The results were analyzed using Fisher Exact test and differences between various groups did not reach statistical significance using two-tailed P values and shown in Table 11 below. The number of African American and Asians analyzed in the study were low to make any definite conclusion regarding frequency of premutation, the lack of detection of premutations in self-identified Asians is consistent with published reports (Otsuka et al. Brain Dev. 2010; 32(2): 110-114; Tzeng et al. Am. J. Med. Genet. A 2005; 133A: 37-43).

Table 11: Statistical analysis of differences in premutation frequencies in various subpopulation

two tailed P values					
	Caucasian	Hispanic	African American/Black	AJ	Asian
Caucasian		>0.2	>0.5	>0.5	>0.2
Hispanic			>0.2	>0.1	>0.2
African American/Black				>0.5	>0.2
AJ					>0.1
Asian					

[0180] To provide for a better approximation of the permutation carrier frequency in US subpopulation or ethnic group, the permutation carrier frequency in each ethnic group in the study was calculated, the data was then extrapolated using the US Census Bureau and approximate Census numbers. Table 12 shows the approximation of permutation carrier frequencies in US subpopulation or ethnic groups based on US Census Bureau data. The permutation carrier frequency in US females is estimated to be 1:177 and that for US males to be 1:396. This calculated frequency may be revised once better information obtained from African Americans and Asians is taken into account.

Table 12: Estimated distribution of permutation carrier frequency in different US subpopulation or ethnic group based on US Census Bureau data

Ethnicity	% of US population ¹	Premutation Carrier Frequency	Estimated Number of females in the U.S. Population ¹	Number of premutation Carriers
Caucasians ²	65.6	1:169	98,759,970 ³	584,379
Hispanic/Latino	15.4	1:287	23,875,899	83,191
Black of African- American	12.8	1:124	19,844,903	160,040
Asian	4.5	1:1500 ³	6,976,724	4,651
Ashkenazi Jewish	2.2	1:134	$2,945,160^3$	21,979
Total			152,402,657	854,239
U.S. Female premutation carrier frequency				1:178
U.S. Male premutation carrier frequency ⁴			=1/178 × (1-0.107) × 0.5	1:400
Full mutation frequency in the U.S. population ⁴			$=1/178 \times 0.107 \times 0.5$	1: 3335

^{1.} Information from U.S. Census Bureau: State and County QuickFacts accessed May 31, 2010 (http://quickfacts.census.gov/qfd/states/00000.html). Numbers may not add to actual total U.S. population due to self identification of two or more ethnicities. Estimates of female numbers assumes females account for 50% of the population.

^{2.} Caucasians did not include Ashkenazi Jews to better estimate overall carrier frequency.

^{3.} Rough estimate from published reports from Taiwan and Japan.

^{4.} Calculated using the equations discussed above.

EXAMPLE 12

Estimation of full mutation frequency in US population

[0181] Full mutation frequency in an individual (male or female) can be calculated based on the premutation carrier frequency of a female and using the following equation:

Full mutation frequency = premutation carrier frequency in females x S x 0.5 where S is the rate of expansion of a premutation into a full mutation and was calculated to be 0.107 (i.e. $\sim 10.7\%$ of premutations expand to full mutations), while 0.5 is the probability of inheriting either allele.

[0182] The calculated full mutation frequency in the U.S. population (males and females) is $1/177 \times 0.107 \times 0.5 = 1:3308$. In Caucasians, the predicted full mutation carrier frequency is 1:3159 individuals while in Hispanics it is 1:5364, and in Ashkenazi Jews 1:2505.

EXAMPLE 13

Estimation of FXTAS and FXPOI incidence in US population

[0183] The incidence of FXTAS and FXPOI in a populations was calculated using the following equation:

incidence of FXTAS = premutation carrier frequency in males x P_{FXTAS} x frequency of alleles ≥ 70 CGG repeats

where P_{EXTAS} is penetrance in patients of larger premutation alleles (\geq 70 CGG repeats) estimated to be 1/3 (0.33).

The incidence of FXPOI incidence was calculated using the following equation:

incidence of FXPOI = premutation carrier frequency in females $x P_{FXPOI} x$ frequency of alleles $\geq 70 CGG$ repeats

where P_{FXPOI} is the penetrance of large premutation alleles(\geq 70 CGG repeats) in FXPOI, estimated to be \sim 20%-28%.

[0184] To determine the incidence of FXTAS in US population, the calculated frequency of premutation in males was first determined based on the premutation carrier frequency in females as discussed above. The calculated premutation carrier frequency in males was 1:396. The calculated incidence of FXTAS in the US-screened population was estimated to be 1:4808 males based on the above equation.

[0185] For incidence of FXPOI, the penetrance of 20% of larger premutations alleles and frequency of alleles ≥70 CGG repeats were used in the calculation similar to FXTAS. Using

this information the calculated incidence of FXPOI in females is 1:3557 in US. This incidence might actually vary if women who suffer from FXPOI were not present in the CF-screened samples.

EXAMPLE 14

Detection of CCG repeats in the 5'-UTR of FMR2 gene

[0186] The CCG repeats in the FMR2 gene can be amplified by PCR. The PCR reaction will include a fluorescently labeled forward primer and reverse primers comprising CGG triplets such that the reverse primers can anneal to CCG repeat sequence. In one example, the reverse primer may comprise the junction sequence of the 3'-end of the CCG repeat tract and flanking sequence and a portion of the sequence downstream of CCG repeat. The reverse primers having at least two CGG triplets and without such junction sequence will hybridize randomly inside the CCG repeat sequence and initiate DNA synthesis. This will result in a range of amplicons of different sizes (stutter). The reverse primers may further comprise a M13 linker sequence at its 5'-end such that a M13 primer can anneal to this linker sequence and further amplify the amplicons thus generated.

[0187] The PCR mix for amplification of all normal, premutation, and full mutations of the CCG repeats may contain 1X PCR buffer; 2mM MgCl₂; 6% DMSO (Sigma); 1.7X Q solution (Qiagen); 0.2 mM each of 7-deaza-2-deoxyGTP, dATP, dCTP, and dTTP; 1 unit of Roche FastStarTaq polymerase; and 0.6 μM each of the forward primer (SEQ ID NO: 76), reverse primer, and M13 reverse linker primer. Several variations of the reverse primers can be tested together with the forward primer. The sequences of the reverse primers are listed in Table 3 above.

[0188] Three microliters of genomic DNA (10-150 ng total) may be added for a final PCR reaction volume of 15 μl. The cycling program, can be performed on an ABI 9700 thermal cycler, can be started by incubating the mix at 98 °C for 10 minutes to activate the polymerase, followed by 10 cycles of denaturation at 97 °C for 35 sec, annealing at 64 °C for 2 minutes, and extension at 68 °C for 8 minutes, followed by 25 cycles of denaturation at 97 °C for 35 sec, annealing at 64 °C for 2 minutes, and extension at 68 °C for 8 min 20 sec (with 20 sec extension each additional cycle). For all 35 cycles, the ramp rate can be adjusted to 64% for the denaturation and extension steps (ramping up at 0.5 °C), and to 25% for the annealing step (ramping down at 0.5 °C). Two microliters of the PCR product (some times diluted 1:5 in H₂O) may be mixed with HiDi Formamide and Map marker-1000 size standard

(Bioventures Inc., Murfreesboro, TN) and samples can be injected on an ABI 3730 DNA Analyzer equipped with a 36-cm capillary loaded with Pop-7 polymer (Applied Biosystems) at 3 volts for 7 seconds. Data can be analyzed using GeneMapper software (Applied Biosystems).

[0189] Unless otherwise defined, all technical and scientific terms used herein have the same meaning as commonly understood by one of ordinary skill in the art to which this invention belongs. All nucleotide sequences provided herein are presented in the 5' to 3' direction.

[0190] The inventions illustratively described herein may suitably be practiced in the absence of any element or elements, limitation or limitations, not specifically disclosed herein. Thus, for example, the terms "comprising", "including," containing", etc. shall be read expansively and without limitation. Additionally, the terms and expressions employed herein have been used as terms of description and not of limitation, and there is no intention in the use of such terms and expressions of excluding any equivalents of the features shown and described or portions thereof, but it is recognized that various modifications are possible within the scope of the invention claimed.

[0191] Thus, it should be understood that although the present invention has been specifically disclosed by preferred embodiments and optional features, modification, improvement and variation of the inventions embodied therein herein disclosed may be resorted to by those skilled in the art, and that such modifications, improvements and variations are considered to be within the scope of this invention. The materials, methods, and examples provided here are representative of preferred embodiments, are exemplary, and are not intended as limitations on the scope of the invention.

[0192] The invention has been described broadly and generically herein. Each of the narrower species and subgeneric groupings falling within the generic disclosure also form part of the invention. This includes the generic description of the invention with a proviso or negative limitation removing any subject matter from the genus, regardless of whether or not the excised material is specifically recited herein.

[0193] In addition, where features or aspects of the invention are described in terms of Markush groups, those skilled in the art will recognize that the invention is also thereby described in terms of any individual member or subgroup of members of the Markush group.

[0194] All publications, patent applications, patents, and other references mentioned herein are expressly incorporated by reference in their entirety, to the same extent as if each were incorporated by reference individually. In case of conflict, the present specification, including definitions, will control.

[0195] Other embodiments are set forth within the following claims.

What is claimed is:

1. A method for determining the presence or absence of an expansion of the CGG repeat tract in the 5'-untranslated region (5'-UTR) of fragile X related mental retardation gene 1 (FMR1) in an individual, said method comprising:

- a) amplifying all or portion of said 5'-UTR CGG repeat tract using as template FMR1 encoding nucleic acid obtained from a sample from said individual, said amplifying being achieved with at least a primer pair, said primer pair comprising
 - (i) a downstream primer that comprises a CCG repeat segment which contains less than four consecutive CCG triplet repeats and hybridizes to at least a portion of said CGG repeat tract and
- (ii) said primer pair comprising an upstream primer that comprises a CGG repeat segment which contains less than four consecutive CGG triplet repeats and hybridizes to at least a portion of the complement of said CGG repeat tract, said amplifying generating amplicons that contain at least a portion of the CGG repeat tract;
- b) separating said amplicons according to size, wherein the size of said amplicons indicates the presence or absence of expansion of said CGG repeat tract in said 5'-UTR of FMR1 gene.
 - 2. The method of claim 1 wherein said expansion of said CGG tract is a full mutation.
 - 3. The method of claim 1 wherein said expansion of said CGG tract is a premutation.
- 4. The method of claim 1 wherein said primer pair comprises a downstream primer that comprises less than four consecutive CCG triplet repeats and hybridizes to at least a portion of said CGG repeat tract and said upstream primer hybridizes upstream of the CGG repeat tract.
- 5. The method of claim 1 wherein said primer pair comprises an upstream primer that comprises less than four consecutive CGG triplet repeats and hybridizes to at least a portion of the complement of said CGG repeat tract, said amplifying generating amplicons that contain at least a portion of the CGG repeat tract and said downstream primer hybridizes downstream of the CGG tract.
- 6. The method of claim 1, wherein said downstream primer comprises a second nucleic acid sequence with no sequence homology to said FMR1 gene, said second nucleic acid

sequence located 5' to said CCG triplet repeat sequence and said upstream primer comprises additional nucleic acid sequence with no sequence homology to said FMR1 gene, said additional nucleic acid sequence located 5' to said CGG repeat sequence.

- 7. The method of claim 6, further comprising a third primer that hybridizes to said second nucleic acid sequence in said downstream primer and/or a fourth primer that hybridizes to said additional nucleic acid sequence in said upstream primer.
- 8. The method of claim 7, wherein said third primer or said fourth primer has a sequence of SEQ ID NO: 11.
- 9. The method of claim 1, wherein said downstream primer further comprises a second nucleic acid sequence 5' to said CCG triplet repeat sequence, said second nucleic acid sequence comprising sequence adjacent to the 3' end of the CGG repeat tract in the FMR1 gene and said upstream primer further comprises a first additional nucleic acid sequence 5' to said CGG repeat sequence, said first additional nucleic acid sequence comprising sequence adjacent to the 5' end of the CGG repeat tract in said FMR1 gene.
- 10. The method of claim 9 wherein said downstream primer hybridizes to the junction between the 3' end of the CGG repeat tract in the FMR1 gene and sequence directly 3' thereto.
- 11. The method of claim 9 wherein said upstream primer hybridizes to the junction between the 5' end of the CGG repeat tract in the FMR1 gene and sequence directly 5' thereto.
- 12. The method of claim 9, wherein said downstream primer further comprises a third nucleic acid sequence 5' to said second sequence, said third sequence having no sequence homology to said FMR1 gene and said upstream primer further comprises a second additional sequence 5' to said first additional sequence, said second additional sequence having no sequence homology to said FMR1 gene.
- 13. The method of claim 12, further comprising a third primer that hybridizes to said third nucleic acid sequence in said downstream primer and/or a fourth primer that hybridizes to said second additional nucleic acid sequence in said upstream primer.

14. The method of claim 13 wherein said third primer or fourth primer has a sequence of SEQ ID NO: 11.

- 15. The method of claim 1, wherein said downstream primer consists essentially of a sequence selected from the group consisting of SEQ ID NOs 3-6, and 14-40, and wherein said upstream primer consists essentially of a sequence selected from the group consisting of SEQ ID NOs: 41-69.
- 16. The method of claim 1, wherein said downstream primer comprises two sets of two consecutive CGG triplets separated from each other by at least one nucleotide and said upstream primer comprises two sets of two consecutive CCG triplets separated from each other by at least one nucleotide.
- 17. The method of claim 16, wherein said downstream primer consists essentially of a sequence selected from the group consisting of SEQ ID NOs: 15-22 and 29-36 or wherein said upstream primer consists essentially of a sequence selected from the group consisting of SEQ ID NOs. 42-49 and 56-66.
- 18. The method of claim 1 wherein said CCG repeat segment in said downstream primer is any of the sequences shown in Figure 15.
- 19. The method of claim 1 wherein said CGG repeat segment in said upstream primer is any of the sequences shown in Figure 16.
- 20. The method of claim 1, wherein said amplification is carried out is a reaction mixture comprising 7-deaza-2'-deoxy GTP.
- 21. The method of claim 20, wherein said 7-deaza-2'-deoxy GTP substitutes for dGTP in said reaction mixture.
- 22. The method of claim 1 wherein said amplification is cyclic and the ramp rate between cycles is selected at about 60-70% (about 0.4°C to 0.6°C/sec) for the denaturation and/or extension steps and about 20-30% (about 0.4°C to 0.6°C/sec) for the annealing step.
 - 23. The method of claim 1, wherein said separation is by capillary electrophoresis.
- 24. The method of claim 1, wherein said separating of amplicons according to size reveals a stutter pattern, said stutter pattern indicating the presence of an expansion of said CGG repeat tract in said 5'-UTR of FMR1 gene.

25. A method for identifying with respect to the CGG tract of the FMR1 gene, a female that is homozygous normal versus a female that is heterozygous and having a normal allele and an allele with expanded CGG tract length, comprising:

- a) amplifying all or portion of said 5'-UTR CGG repeat tract using as template FMR1 encoding nucleic acid obtained from a sample from said female, said amplifying being achieved with at least a primer pair, said primer pair comprising a downstream primer that comprises a CCG repeat segment which contains less than four consecutive CCG triplet repeats and hybridizes to at least a portion of said CGG repeat tract and/or said primer pair comprising an upstream primer that comprises a CGG repeat segment which contains less than four consecutive CGG triplet repeats and hybridizes to at least a portion of the complement of said CGG repeat tract, said amplifying generating amplicons that contain at least a portion of the CGG repeat tract;
- b) separating said amplicons according to size so as to estimate the number of CGG repeats in said amplicons, wherein detection of only a normal number of CGG repeats indicates a homozygous normal female and detection of both amplicons with a normal number of CGG repeats and detections of amplicons with an abnormally high number of CGG repeats indicates a female heterozygous and having an allele with an expanded CGG tract length.
- 26. The method of claim 25, wherein a normal number of CGG repeats is between 4 to 44 repeats while an abnormally high number of repeats is greater than 54 repeats.
- 27. The method of claim 25 wherein said expanded CGG tract length represents a premutation with CGG repeats of between 55-200
- 28. The method of claim 25 wherein said expanded CGG tract length represents a full mutation with CGG tract represents of greater than 200.
- 29. A method of screening newborns to determine if they have an expanded CGG tract in the FMR1 gene, comprising
- a) amplifying all or portion of said 5'-UTR CGG repeat tract using as template FMR1 encoding nucleic acid obtained from a sample from said newborn, said amplifying being achieved with at least a primer pair, said primer pair comprising a downstream primer that comprises a CCG repeat segment which contains less than four consecutive CCG triplet repeats and hybridizes to at least a portion of said CGG repeat tract and/or said primer pair

comprising an upstream primer that comprises a CGG repeat segment which contains less than four consecutive CGG triplet repeats and hybridizes to at least a portion of the complement of said CGG repeat tract, said amplifying generating amplicons that contain at least a portion of the CGG repeat tract;

- b) separating said amplicons according to size so as to estimate the number of CGG repeats in said amplicons, wherein detection of only a normal number of CGG repeats indicates a newborn that has normal CGG tract length FMR1 gene alleles and detection of amplicons with an abnormally high number of CGG repeats indicates a newborn with an allele with an expanded CGG tract length.
- 30. The method of claim 29, wherein a normal number of CGG repeats is between 4 to 44 repeats while an abnormally high number of repeats is greater than 54 repeats.
- 31. The method of claim 29 wherein said expanded CGG tract length represents a premutation with CGG repeats of between 55-200
- 32. The method of claim 29 wherein said expanded CGG tract length represents a full mutation with CGG tract represents of greater than 200.
- 33. A method of estimating premutation carrier frequency in FMR gene in males in a population or ethnic group, said method comprising:
- a) determining the total number of premutation in FMR gene in females representative of the population or ethnic group using the method of claim 1;
- b) calculating the premutation carrier frequency in FMR gene in females in said population or ethnic group by dividing the total number of females identified as having said premutation by the total number of individuals tested in said population or ethnic group;
- c) estimating said premutation carrier frequency in FMR gene in males in said population or ethnic group using the formula: premutation carrier frequency in males = premutation carrier frequency in females X (1-S) X 0.5, wherein S is the rate of expansion of a premutation into a full mutation in said population or ethnic group.
- 34. The method of claim 33, wherein said rate of expansion of a premutation into a full mutation in US population is 0.107.

35. A method of estimating the full mutation frequency in FMR gene in a population or ethnic group, said method comprising:

- a) determining the total number of premutation in FMR gene in females representative of said population or ethnic group using the method of claim 1;
- b) calculating the premutation carrier frequency in FMR gene in females in said population or ethnic group by dividing the total number of females identified as having said premutation by the total number of individuals tested in said population or ethnic group;
- c) estimating said full mutation frequency in FMR gene in said population or ethnic group using the formula: full mutation carrier frequency = premutation carrier frequency in females X S X 0.5, wherein S is the rate of expansion of a premutation into a full mutation in said population or ethnic group.
- 36. The method of claim 35, wherein said rate of expansion of a premutation into a full mutation in US population is 0.107.
- 37. A method of estimating the prevalence of Fragile X-associated tremor and ataxia syndrome (FXTAS) in males in a population or ethnic group, said method comprising:
- a) determining the total number of premutation in FMR gene in females representative of the population or ethnic group using the method of claim 1;
- b) calculating the premutation carrier frequency in FMR gene in females in said population or ethnic group by dividing the total number of females identified as having said premutation by the total number of individuals tested in said population or ethnic group;
- c) estimating said premutation carrier frequency in FMR gene in males in said population or ethnic group using the formula: premutation carrier frequency in males = premutation carrier frequency in females X (1-S) X 0.5, wherein S is the rate of expansion of a premutation into a full mutation in said population or ethnic group;
- d) estimating the prevalence of FXTAS in males in said population or ethnic group using the formula: prevalence of FXTAS = premutation carrier frequency in males X $P_{\rm FXTAS}$ X frequency of alleles ≥ 70 CGG repeats, wherein $P_{\rm FXTAS}$ is the penetrance in patients having permutation ≥ 70 CGG repeats.

38. The method of claim 37, wherein said penetrance in male patients having permutation \geq 70 repeats in US population is 0.33.

- 39. A method of estimating the prevalence of Fragile X-associated primary ovarian insufficiency (FXPOI) in females in a population or ethnic group, said method comprising:
- a) determining the total number of premutation in FMR gene in females representative of the population or ethnic group using the method of claim 1;
- b) calculating the premutation carrier frequency in FMR gene in females in said population or ethnic group by dividing the total number of females identified as having said premutation by the total number of individuals tested in said population or ethnic group;
- c) estimating the prevalence of FXPOI in females in said population or ethnic group using the formula: prevalence of FXPOI = premutation carrier frequency in females X P_{FXPOI} X frequency of alleles \geq 70 CGG repeats, wherein P_{FXPOI} is the penetrance in patients having permutation \geq 70 CGG repeats.
- 40. The method of claim 39, wherein said penetrance in female patients having permutation ≥70 repeats in US population is 0.2.
 - 41. The method of claim 20, wherein said reaction mixture further comprises GTP.
- 42. The method of claim 41, wherein the concentration of GTP in said reaction mixture is 0.2 mM.
- 43. A method for determining the presence or absence of an expansion of the CCG repeat tract in the 5'-untranslated region (5'-UTR) of fragile X related mental retardation gene 2 (FMR2) in an individual, said method comprising:
- a) amplifying all or portion of said 5'-UTR CCG repeat tract using as template FMR2 encoding nucleic acid obtained from a sample from said individual, said amplifying being achieved with at least a primer pair, said primer pair comprising a downstream primer that comprises one or more CCG repeat segment and hybridizes to at least a portion of said CCG repeat tract and/or said primer pair comprising an upstream primer that comprises one or more CCG triplets and hybridizes to at least a portion of the complement of said CCG repeat tract, said amplifying generating amplicons that contain at least a portion of the CCG repeat tract;

b) separating said amplicons according to size, wherein the size of said amplicons indicates the presence or absence of expansion of said CCG repeat tract in said 5'-UTR of FMR2 gene.

- 44. The method of claim 43 wherein said expansion of said CGG tract is full mutation.
- 45. The method of claim 43 wherein said expansion of said CGG tract is premutation.

(1/25)

Figure 1 (SEQ ID NO:1)

```
13081 aatccatqtc ccttaaaqqq cacaqqqtqt ctccacaqqq ccqcccaaaa tctqqtqaqa
13141 gagggcgtag acqcctcacc ttctgcctct acgggtcaca aaagcctggg tcaccctggt
13201 tgccactgtt cctagttcaa agtcttcttc tgtctaatcc ttcaccccta ttctcgcctt
13321 ccaaaccaaa ccaaaccaga ccagacaccc cctcccgcgg aatcccagag aggccgaact
13381 gggataaccg gatgcatttg atttcccacg ccactgagtg cacctctgca gaaatgggcg
13441 ttctggccct cgcgaggcag tgcgacctgt caccgccctt cagccttccc gccctccacc
13501 aagcccgcgc acgcccggcc cgcgcgtctg tctttcgacc cggcaccccg gccggttccc
13561 agcagegege atgegegege teccaggeca ettgaagaga gagggegggg eegagggget
13621 gagcccgcgg ggggagggaa cagcgttgat cacgtgacgt ggtttcagtg tttacacccg
13681 cagcgggccg ggggttcggc ctcagtcagg cgctcagctc cgtttcggtt tcacttccgg
13741 tggagggccg cctctgagcg ggcggcggc cgacggcgag cgcgggcggc ggcggtgacg
13801 gaggegeege tgecaggggg egtgeggeag egeggeggeg geggeggeg eggeggegge
13861 ggaggeggeg geggeggeg eggeggegge ggetgggeet egagegeeeg eageceacet
13921 ctcgggggcg ggctcccggc gctagcaggg ctgaagagaa gatggaggag ctggtggtgg
13981 aagtgegggg etecaatgge getttetaea aggtaettgg etetagggea ggeeceatet
14041 tegecettee tteecteect tttettettg gtgteggegg gaggeaggee eggggeecte
14101 ttcccgagca ccgcgcctgg gtgccagggc acgctcggcg ggatgttgtt gggagggaag
14161 gactggactt ggggcctgtt ggaagcccct ctccgactcc gagaggccct agcgcctatc
14221 gaaatgagag accagcgagg agagggttct ctttcggcgc cgagccccgc cggggtgagc
14281 tggggatggg cgagggccgg cggcaggtac tagagccggg cgggaagggc cgaaatcggc
14341 gctaagtgac ggcgatggct tattccccct ttcctaaaca tcatctccca gcgggatccg
14401 ggcctgtcgt gtgggtagtt gtggaggagc gggggggcgct tcagccgggc cgcctcctgc
14461 agegecaaga gggetteagg teteetttgg ettetetttt eeggtetage attgggaett
14521 cggagagctc cactgttctg ggcgagggct gtgaagaaag agtagtaaga agcggtagtc
14581 ggcaccaaat cacaatggca actgattttt agtggcttct ctttgtggat ttcggaggag
14641 attttagatc caaaagtttc aggaagaccc taacatggcc cagcagtgca ttgaagaagt
14701 tgatcatcgt gaatattcgc gtcccccttt ttgttaaacg gggtaaattc aggaatgcac
14761 atgcttcagc gtctaaaacc attagcagcg ctgctactta aaaattgtgt gtgtgttt
14821 aagtttccaa agacctaaat atatgccatg aaacttcagg taattaactg agagtatatt
14881 attactaggg cattttttt ttaactgagc gaaaatattt ttgtgcccct aagaacttga
14941 ccacatttcc tttgaatttg tggtgttgca gtggactgaa ttgttgaggc tttatatagg
15001 cattcatggg tttactgtgc tttttaaagt tacaccattg cagatcaact aacacctttc
15061 agttttaaaa ggaagattta caaatttgat gtagcagtag tgcgtttgtt ggtatgtagg
15121 tgctgtataa attcatctat aaattctcat ttccttttga atgtctataa cctctttcaa
15181 taatatccca ccttactaca gtattttggc aatagaaggt gcgtgtggaa ggaaggctgg
15241 aaaatagcta ttagcagtgt ccaacacaat tcttaaatgt attgtagaat ggcttgaatg
15301 tttcagacag gacacgtttg gctataggaa aataaacaat tgactttatt ctgtgtttac
15361 caattttatq aaqacatttq qaqatcaqta tatttcataa atqaqtaaaq tatqtaaact
15421 gttccatact ttgagcacaa agataaagcc ttttgctgta aaaggaggca aaaggtaacc
15481 ccgcgtttat gttcttaaca gtctcatgaa tatgaaattg tttcagttga ctctgcagtc
15541 aaaattttaa tttcattgat tttattgatc cataatttct tctggtgagt ttgcgtagaa
15601 tcgttcacgg tcctagatta gtggttttgg tcactagatt tctggcacta ataactataa
15661 tacatataca tatatatgtg tgagtaacgg ctaatggtta ggcaagattt tgattgacct
15721 gtgatataaa cttagattgg atgccactaa agtttgctta tcacagaggg caagtagcac
15781 attatggcct tgaagtactt attgttctct tccagcaact tatgatttgc tccagtgatt
15841 ttgcttgcac actgactgga atataagaaa tgccttctat ttttgctatt aattccctcc
15901 ttttttgttt tgttttgtaa cgaagttgtt taacttgaag gtgaatgaag aataggttgg
15961 ttgcccctta gttccctgag gagaaatgtt aatacttgaa caagtgtgtg tcagacaaat
16021 tgctgttatg tttatttaat taagtttgat ttctaagaaa atctcaaatg gtctgcactg
16081 atggaagaac agtttctgta acaaaaaagc ttgaaatttt tatatgactt ataatactgc
16141 tgtgagtttt aaaagtaaag caaaagtaaa ctgagttgct tgtccagtgg gatggacagg
16201 aaagatgtga aataaaaacc aatgaaaaat gaactgctgt ggagaagtgt tacatttatg
16261 gaaaaagaaa taggaacctt gttcatcaaa ttgatagaaa agcttttaaa actaaacaaa
```

(2/25)

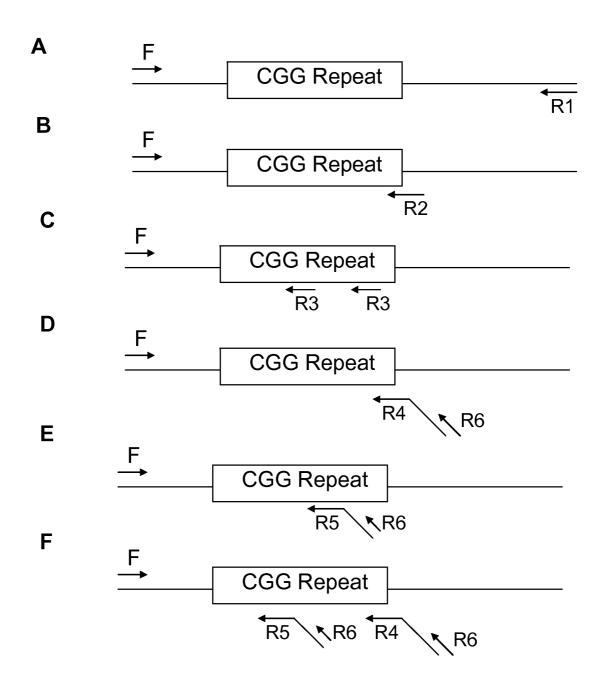


FIGURE 2 (page 1 of 2)

(3/25)

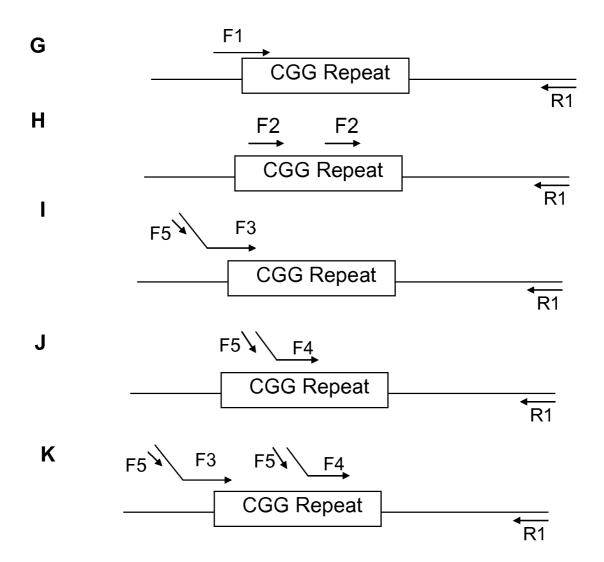


FIGURE 2 CONTINUED (page 2 of 2)

(4/25)

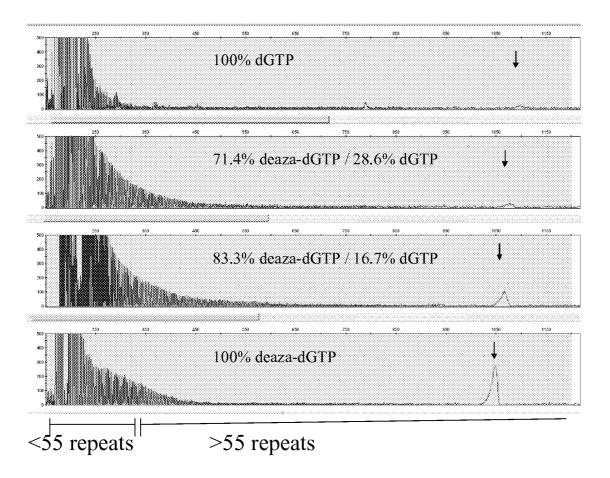


FIGURE 3

(5/25)

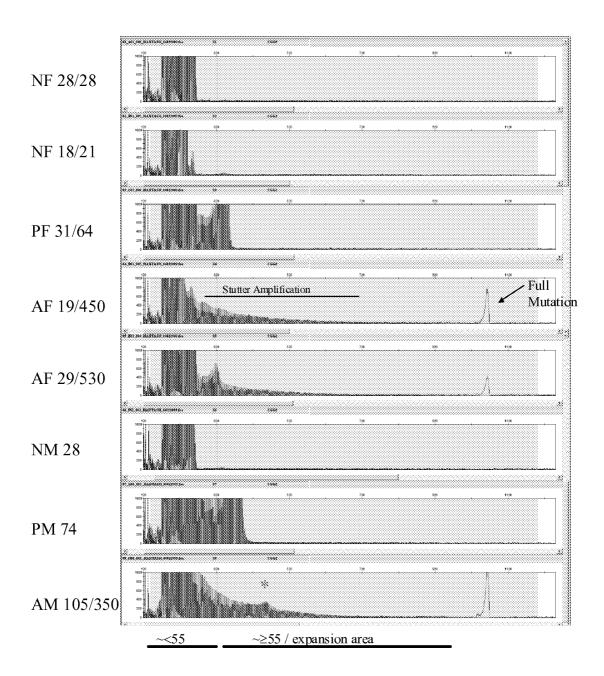


FIGURE 4

(6/25)

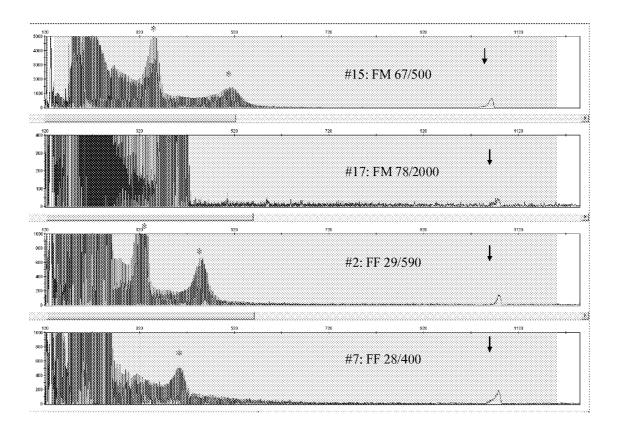


FIGURE 5

(7/25)

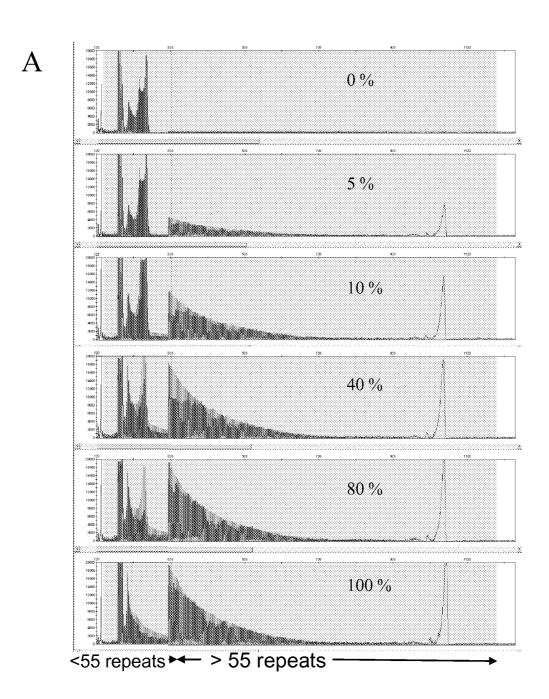


FIGURE 6 (Page 1 of 2)

(8/25)

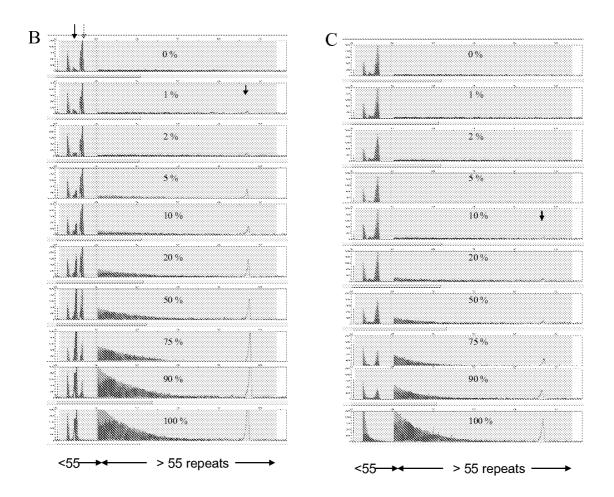


FIGURE 6 CONTINUED (Page 2 of 2)

(9/25)


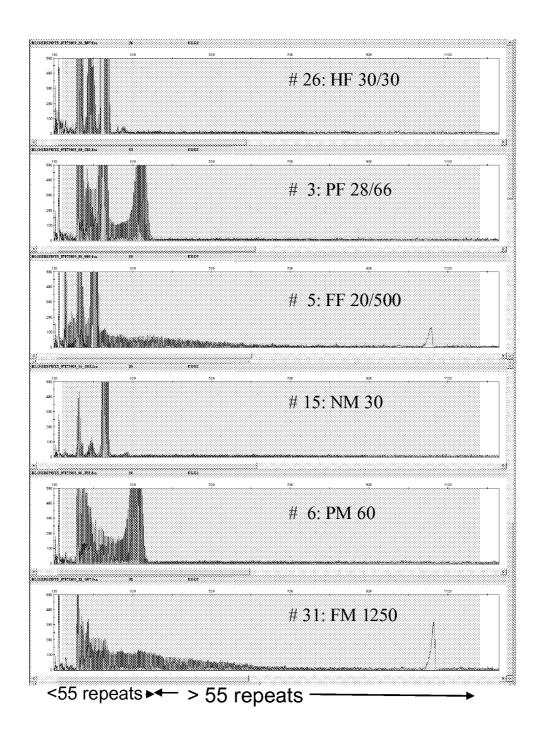


FIGURE 7


(10/25)

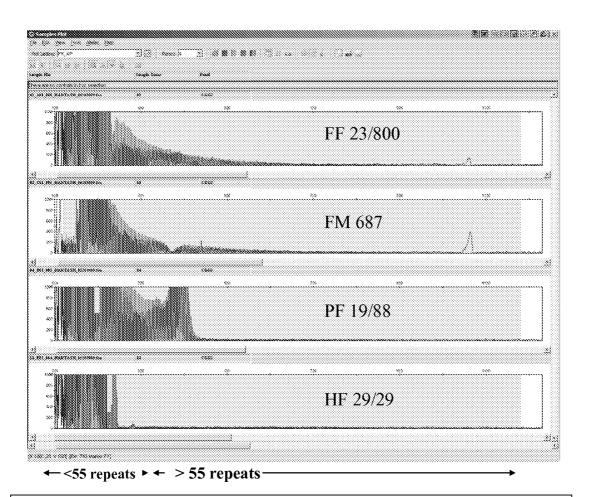

FMR1R: CAGGAAACAGCTATGACC ctcgaggcccagCCGCCGCCCC

FIGURE 8


(11/25)

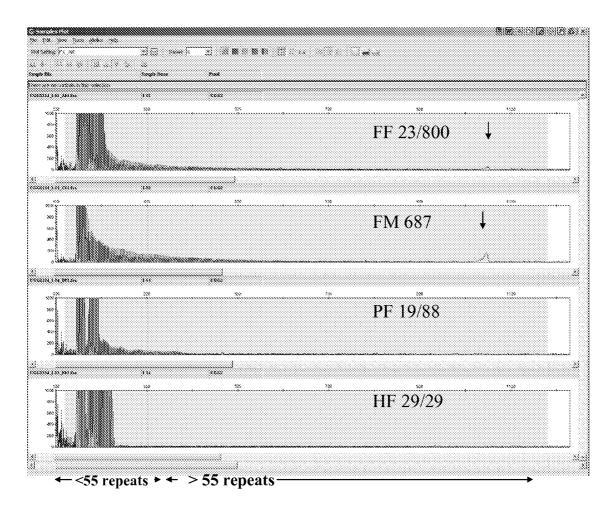

CGG1: CAGGAAACAGCTATGACC CG CCG CCG CC

FIGURE 9


(12/25)

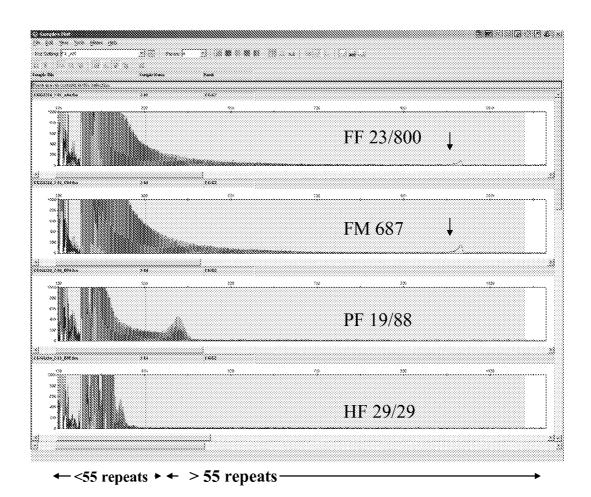

CGG2: CAGGAAACAGCTATGACC CCG CCG XCG CCG CCG CC

FIGURE 10


(13/25)

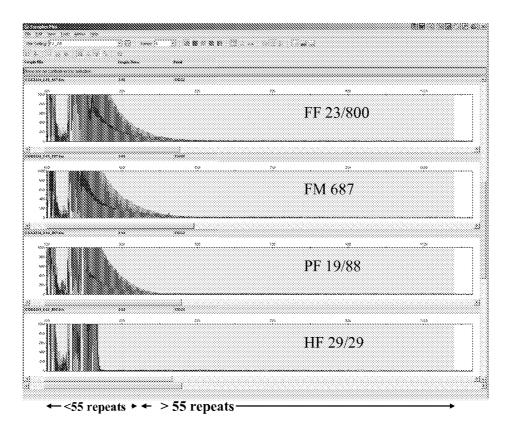

CGG3: CAGGAAACAGCTATGACCCCGCCGCCGCCGCCGCCGCCGCCGCCGCCG

FIGURE 11

(14/25)

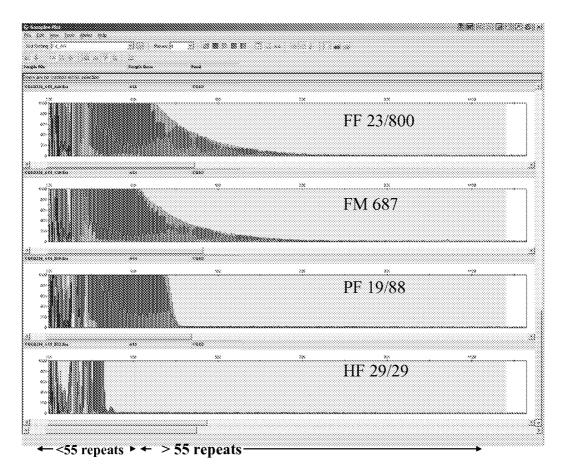

CGG4: CAGGAAACAGCTATGACC CCG CCG CCG CCG CCG CCG

FIGURE 12

(15/25)


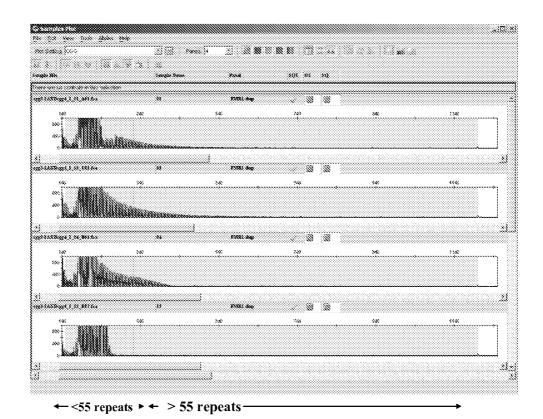


FIGURE 13

(16/25)

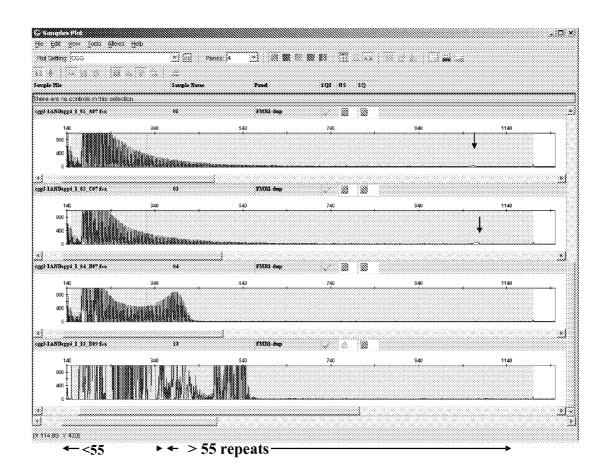

CGG4-1: CAGGAAACAGCTATGACC CCG CCG CCG CCG CCG

FIGURE 14

(17/25)

Primer Sequence	Variable Bases and Lengths	SEQ ID NO:
(CCG)n	n= 1-15	78
(CCG) ₂ DCG (CCG) ₂ CC	D= A, G, T	79
CG (CCG) ₂ DCG (CCG) ₂ CC	D= A, G, T	80
(CCG) ₂ DCG (CCG) ₂	D= A, G, T	81
CG(CCG) ₂ DCG (CCG) ₂ C	D= A, G, T	82
G(CGG) ₂ DGG (CGG) ₂ C	D= A, G, T	83
(CCG) ₂ N (CCG) ₂	N= A, C, G, T	84
CG(CCG) ₂ N (CCG) ₂ CC	N= A, C, G, T	85
G(CCG) ₂ N (CCG) ₂ CC	N= A, C, G, T	86
CG(CCG) ₂ N (CCG) ₂ C	N= A, C, G, T	87
(CCG) ₂ DCG (CCG) ₂ DCG (CCG) ₂ CC	D= A, G, T	88
CG(CCG) ₂ DCG (CCG) ₂ DCG (CCG) ₂ CC	D= A, G, T	89
(CCG) ₂ N (CCG) ₂ N (CCG) ₂	N= A, C, G, T	90
(CCG) ₂ N (CCG) ₂ N (CCG) ₂ CC	N= A, C, G, T	91
CG (CCG) ₂ N (CCG) ₂ N (CCG) ₂ CC	N= A, C, G, T	92
CGCCGNCCGCC	N= A, C, G, T	93
GCCGNCCGCC	N= A, C, G, T	94
CGCCGCCGCC		95

FIGURE 15

(18/25)

Primer Sequence	Variable Bases and Lengths	SEQ ID NO:
(CGG)n	n= 1-15	96
(CGG) ₂ DGG (CGG) ₂ CG	D= A, G, T	97
GG (CGG) ₂ DGG (CGG) ₂ CG	D= A, G, T	98
(CGG) ₂ DGG (CGG) ₂	D= A, G, T	99
GG(CGG) ₂ DGG (CGG) ₂ C	D= A, G, T	100
G(CGG) ₂ DGG (CGG) ₂ C	D= A, G, T	101
(CGG) ₂ N (CGG) ₂	N= A, C, G, T	102
GG(CGG) ₂ N (CGG) ₂ CG	N= A, C, G, T	103
G(CGG) ₂ N (CGG) ₂ CG	N= A, C, G, T	104
G(CGG) ₂ N (CGG) ₂ C	N= A, C, G, T	105
(CGG) ₂ DGG (CGG) ₂ DGG (CGG) ₂ CG	D= A, G, T	106
GG(CGG) ₂ DGG (CGG) ₂ DGG (CGG) ₂ CG	D= A, G, T	107
(CGG) ₂ N (CGG) ₂ N (CGG) ₂	N= A, C, G, T	108
(CGG) ₂ N (CGG) ₂ N (CGG) ₂ CG	N= A, C, G, T	109
GG (CGG) ₂ N (CGG) ₂ N (CGG) ₂ CG	N= A, C, G, T	110
GGCGGNCGGCG	N= A, C, G, T	111
GCGGNCGGCG	N= A, C, G, T	112
GGCGGCGCG	N= A, C, G, T	113

FIGURE 16

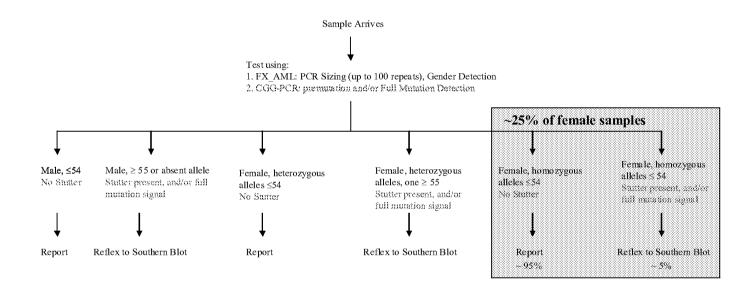
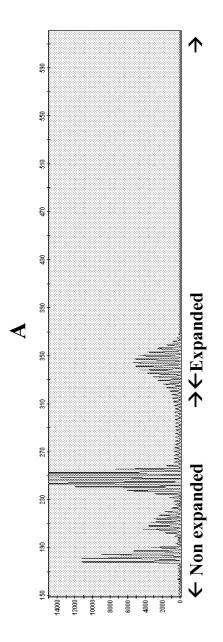


FIGURE 17




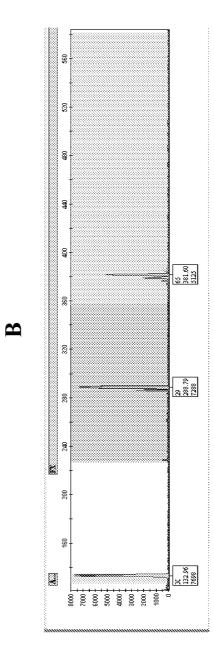


FIGURE 18

WO 2011/097503 (ZI/ZJ) PCT/US2011/023777

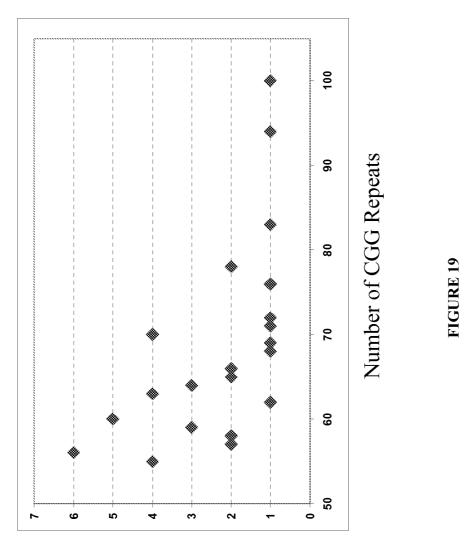

Frequency of CGG Repeats

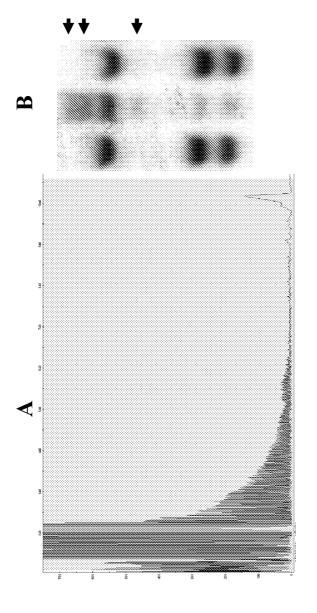



FIGURE 20

(23/25)

75) (SEQ ID NO: 61 egeogetgee geceeggetg eegegeegeg eegetgeete tgeeeeeggee geeeeegeeg 181 geggetgeeg eegeagegte gggtegetgg gtgegeggge taeegeggae egageggaee 241 cgagtgggcg accaggcgct tgcccgccca gtgccactgc cgccgcttcc tcgccggagc 121 ccgctgccgc cgccggcccg cagccagcca ggcgggcggc ccagcccgcc tgagcccgca 301 acaggaccag acacetecag egecegetge tgetgeegat geggeeegga eaettttage 361 tgggcgggag ggctggagag ccgggggccg ccgagaaccg ccagcgagct gtgccgagag 421 ccgcgccgac ccgctgcgat cagggacagg cgcccgccg ccgccgccgc ctggccgct

FIGURE 21

(24/25)

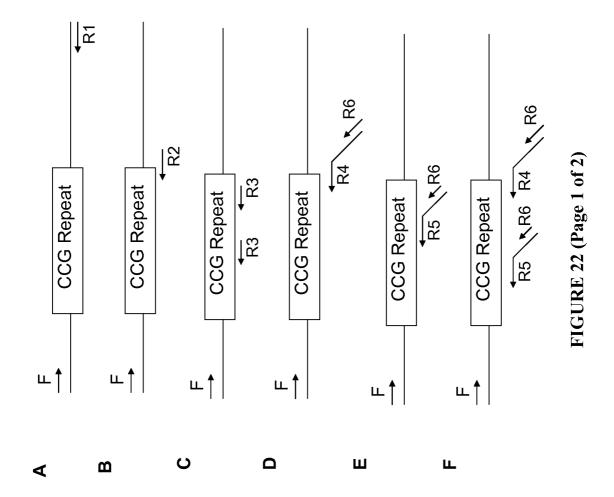

(25/25)

